ULTIMATE 3D PHOTON LED MASK

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PALSAR7 - ABOUT THE BRAND

PALSAR7 IS A BRAND SPECIALIZING IN INNOVATIVE PRODUCTS FOR PROFESSIONAL SKIN CARE FROM THE COMFORT OF YOUR HOME. IT OFFERS MODERN TECHNOLOGIES THAT SUPPORT SKIN HEALTH AND BEAUTY WITHOUT THE NEED TO VISIT A BEAUTY SALON. THE BRAND EMPHASIZES QUALITY, SAFETY AND EASE OF USE OF ITS PRODUCTS. THANKS TO THIS, YOU CAN INDULGE IN PROFESSIONAL CARE ANYTIME, ANYWHERE.

ABOUT THE PRODUCT

THE DEVICE USES THE PRINCIPLE OF LIGHT THERAPY TO ACTIVATE SKIN CELLS AND IMPROVE SKIN CONDITION THROUGH SPECIALLY ADJUSTED LIGHT OF SPECIFIC WAVELENGTHS. IT IS A NON-INVASIVE, PAINLESS AND SAFE METHOD OF SKIN CARE WITHOUT THE USE OF CHEMICALS. AFTER BEING ABSORBED BY THE SKIN, THE LIGHT IS CONVERTED INTO INTRACELLULAR ENERGY, WHICH TRIGGERS A PHOTOCHEMICAL AND ENZYMATIC REACTION. DURING THIS PROCESS, A LARGE AMOUNT OF GLYCOGEN AND PROTEIN IS PRODUCED, WHICH PROMOTE CELL METABOLISM AND COLLAGEN REGENERATION. THANKS TO THIS, EFFECTS SUCH AS REJUVENATION, ACNE TREATMENT, DARK SPOT REDUCTION, SKIN FIRMING, WRINKLE SMOOTHING AND SKIN WHITENING, CAN BE ACHIEVED.

PRODUCT DESCRIPTION

PRODUCT FUNCTIONS

THIS LED MASK USES 99 PIECES OF LEDS IN 4 COLORS (RED, INFRARED, BLUE, YELLOW), THANKS TO THE COMBINATION OF LIGHTS, THE MASK OFFERS 5 MODES. THEIR FUNCTIONS ARE AS FOLLOWS:

RED LIGHT (630 NM \pm 10 NM): THIS LIGHT, CALLED BIOACTIVE LIGHT, IS CHARACTERIZED BY HIGH PURITY, INTENSITY AND UNIFORM DENSITY, WHICH MAKES IT SUITABLE FOR SKIN CARE AND HEALTH. RED LIGHT INCREASES CELL VITALITY, ACCELERATES BLOOD CIRCULATION, PROMOTES CELL METABOLISM, STIMULATES THE REGENERATION OF COLLAGEN AND FIBROUS TISSUE, THEREBY IMPROVING SKIN ELASTICITY AND REDUCING DULLNESS. IT IS A GREAT ANTI-AGING AID, ACTS AS AN ANTIOXIDANT AND PROMOTES SKIN RENEWAL.

BLUE LIGHT (450 NM \pm 10 NM): THIS LIGHT CORRESPONDS TO THE ABSORPTION MAXIMUM OF PORPHYRIN, A METABOLITE OF PROPIONIBACTERIUM ACNES. DURING THE DE-EXCITATION PROCESS AFTER STIMULATION OF PORPHYRIN, A LARGE AMOUNT OF SINGLET ACTIVE OXYGEN IS GENERATED, WHICH CREATES A HIGHLY OXIDATIVE ENVIRONMENT THAT DESTROYS PROPIONIBACTERIUM ACNES. THIS HELPS IN THE TREATMENT OF ACNE.

INFRARED LIGHT/INVISIBLE LIGHT (850 NM \pm 10 NM): THIS LIGHT PROMOTES COLLAGEN PRODUCTION IN THE SKIN, SPEEDS UP ITS METABOLISM, BRIGHTENS SKIN TONE, TIGHTENS PORES, AND HELPS ACHIEVE A MORE YOUTHFUL APPEARANCE. SINCE INFRARED LIGHT IS INVISIBLE, THE EMITTED LIGHT CANNOT BE SEEN WHEN YOU SELECT THIS MODE.

YELLOW LIGHT (580 nm \pm 10 nm): THIS LIGHT HELPS TO IMPROVES THE OXYGEN EXCHANGE FUNCTION OF SKIN CELLS, PROMOTES MICROCIRCULATION, BREAKS DOWN AND LIGHTENS PIGMENT, STIMULATES CELL GROWTH, AND THUS SMOOTHES THE SKIN, BRIGHTENS SKIN TONE, REGENERATES THE SKIN, AND REDUCES WRINKLES.

PURPLE LIGHT (COMBINATION OF RED AND BLUE LIGHTS) (420 nm \pm 10 nm): THIS LIGHT IS A GOOD HELPER IN TREATING ACNE AND REMOVING ACNE SCARS.

HOW TO STORE

- AFTER USE, WIPE THE DEVICE WITH A DRY CLOTH. DO NOT RINSE UNDER WATER. DO NOT SOAK.
- STORE THE DEVICE IN A COOL, SHADY PLACE AND AVOID DIRECT SUNLIGHT OR HIGH TEMPERATURES.
- KEEP THE DEVICE OUT OF THE REACH OF CHILDREN TO PREVENT DAMAGE OR INJURY.
- NEVER ATTEMPT TO OPEN OR REPAIR THE DEVICE YOURSELF. IN THE EVENT OF A MALFUNCTION, CONTACT THE CUSTOMER SERVICE CENTER OR THE SELLER.
- AFTER USE, STORE THE DEVICE IN THE BOX.
- DO NOT PLACE ANY WEIGHT ON THE DEVICE. DO NOT PLACE IT UNDER OTHER OBJECTS AND BE CAREFUL TO AVOID MECHANICAL DAMAGE.
- WE RECOMMEND STORING IT IN A HORIZONTAL POSITION, INCLUDING THE CABLE. NEVER TRY TO DISCONNECT, DEFORM OR TWIST THE CABLE ATTACHED TO THE MASK.

CHARGING

- **PLEASE FULLY CHARGE THE DEVICE BEFORE USING IT FOR THE FIRST TIME.**
- CONNECT THE CONTROLLER TO THE APPROPRIATE ADAPTER (NOT INCLUDED) USING THE USB-C CABLE THAT WAS INCLUDED IN THE PACKAGE.
- **USE ONLY THE BOTTOM PORT FOR CHARGING!**
- PLUG INTO THE POWER OUTLET. WHILE CHARGING, THE SCREEN WILL LIGHT UP AND A BATTERY SYMBOL WILL APPEAR WITH THE CORRESPONDING CHARGE LEVEL.
- A FULL BATTERY SYMBOL INDICATES A FULL CHARGE. ONCE CHARGED, UNPLUG THE ADAPTER FROM THE POWER OUTLET AND YOU CAN START USING THE CONTROLLER AFTER CONNECTING IT TO THE MASK.

HOW TO USE

1. WASH YOUR FACE AND GENTLY PAT DRY WITH A SOFT TOWEL. THEN APPLY TONER TO YOUR FACE.


2. THEN YOU CAN APPLY YOUR FAVORITE NATURAL COSMETICS (SERUM, MASK OR ESSENCE).


3. TAKE OUT THE MASK. CONNECT THE USB-C CABLE CONNECTOR ON THE MASK TO THE PORT ON THE CONTROLLER. USE THE TOP PORT ONLY!


4. HOLD THE MASK CLOSE TO YOUR FACE. GRASP BOTH ENDS OF THE ELASTIC STRAPS AND FASTEN WITH THE VELCRO AT THE BACK OF YOUR HEAD. ADJUST THE ELASTIC STRAPS TO FIT SNUGLY AND PREVENT LIGHT FROM ENTERING THROUGH THE SILICONE EYE PROTECTION.

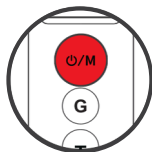5. PRESS AND HOLD THE POWER BUTTON FOR 1.5 SECONDS TO TURN ON THE DEVICE. AFTER POWERING ON, THE DEFAULT MODE IS RED LIGHT. SHORT PRESS THE MODE BUTTON TO SWITCH BETWEEN MODES. THERE ARE FIVE MODES AVAILABLE: 1 - RED LIGHT (R), 2 - BLUE LIGHT (B), 3 - YELLOW LIGHT (Y), 4 - INFRARED LIGHT (IR), 5 - PURPLE LIGHT (R+B). YOU CAN ALSO SWITCH BETWEEN MODES DURING USE BY PRESSING THE MODE BUTTON.

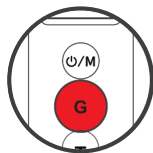6. SHORT PRESS TO ADJUST THE LIGHT INTENSITY. THERE ARE THREE BRIGHTNESS LEVELS: 50%, 75%, 100%.

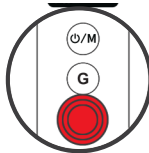7. SHORT PRESS THE TIME SETTING BUTTON TO SELECT THE PROCEDURE TIME. (5/10/15/20/25/30 MIN.)

8. AFTER YOU ARE DONE USING THE MASK, APPLY CREAM TO YOUR FACE TO PREVENT DRYNESS.

RECOMMENDATION:

THE RECOMMENDED USAGE TIME IS 10-30 MINUTES 3-4 TIMES A WEEK. USE THE BLUE LIGHT FOR A MAXIMUM OF 10 MINUTES PER TREATMENT. DO NOT EXCEED THE MAXIMUM WEEKLY USAGE TIME (2 HOURS).

PRODUCT FEATURES

EFFECTIVE

99 LEDS PROVIDE LARGE-AREA HIGH-ENERGY PHOTON SKIN CARE TREATMENT.

VERSATILE

5 LIGHT BEAUTY MODES TO ADDRESS VARIOUS SKIN CONCERNS.

SAFE

NON-INVASIVE SKIN CARE METHOD ENSURES CHEMICAL-FREE AND PAINLESS TREATMENT.

EASY TO USE

ONE-STEP THERAPY AND A BUILT-IN TIMER MAKE THIS MASK EASY TO USE.

COMFORTABLE

SOFT AND COMFORTABLE SILICONE MATERIAL ENSURES PLEASANT USER EXPERIENCE.

FLEXIBLE

THE MASK IS ONLY ATTACHED TO THE CONTROLLER, ALLOWING YOU TO MOVE FREELY DURING USE.

USABLE WITH COSMETICS

USE WITH COSMETIC PRODUCTS SUCH AS ESSENCE, SERUM OR MASK TO IMPROVE EFFECTIVENESS.

TIMING

6 TIMING MODES (5 MIN., 10 MIN., 15 MIN., 20 MIN., 25 MIN. AND 30 MIN.).

INTENSITY SETTINGS

3 ADJUSTABLE LIGHT INTENSITY LEVELS (50%, 75% AND 100%).

CLEAN

WIPEABLE, SILICONE DOES NOT TRANSFER BACTERIA, IS HYPOALLERGENIC AND EASY TO MAINTAIN.

WARNINGS

- THE PRODUCT IS NOT SUITABLE FOR CHILDREN.
- DO NOT USE THE DEVICE IF YOU KNOW THAT YOU ARE PARTICULARLY SENSITIVE TO LIGHT AND HAVE NOT RECENTLY UNDERGONE A LIGHT SENSITIVITY TEST.
- DO NOT USE THE DEVICE IF YOU ARE PREGNANT AS THE RISKS ARE UNKNOWN AND HAVE NOT BEEN ESTABLISHED, CONSULT YOUR DOCTOR.
- DO NOT USE THE DEVICE IF YOU SUFFER FROM EPILEPSY AS SOME LIGHTS OR FREQUENCIES MAY TRIGGER A SEIZURE.
- THE DEVICE IS NOT WATERPROOF. DO NOT ALLOW WATER OR OTHER LIQUIDS TO ENTER ANY PART OF THE DEVICE. THIS MAY DAMAGE IT.
- DO NOT USE THE DEVICE IF YOU ARE TAKING MEDICATION THAT MAY MAKE YOUR SKIN SENSITIVE TO LIGHT.
- DO NOT USE THE DEVICE WHILE SLEEPING, WALKING OR DRIVING.
- USE THE DEVICE ACCORDING TO THE CORRECT INSTRUCTIONS AND AVOID EXCESSIVE USE.
- DISCONTINUE USE IF YOU EXPERIENCE DISCOMFORT OR IF REDNESS OR SKIN DISCOLORATION PERSISTS FOR MORE THAN 24 HOURS. CONSULT A DERMATOLOGIST.
- DO NOT TOUCH THE LEDs WHEN THE DEVICE IS TURNED ON.
- DO NOT USE THE DEVICE IF IT IS VISIBLY DAMAGED.
- DO NOT SUBJECT THE DEVICE TO SHOCKS TO AVOID DEFORMATION OR DAMAGE.
- DO NOT USE CHEMICALS, WATER, OR OTHER CLEANING OR LIQUID PRODUCTS FOR CLEANING.
- DO NOT EXCEED THE RECOMMENDED USAGE TIME.

SPECIFICATIONS

RATED CURRENT: C5V == 1A

RATED VOLTAGE: DC5V

CONTROLLER SIZE: 45X35X135 MM

CONTROLLER WEIGHT: 105 G

BATTERY CAPACITY: 2000MAH

MASK SIZE: 225X240X135 MM

MASK THICKNESS WHEN IN USE: 6 MM

MASK WEIGHT: 380 G

PACKAGING CONTAINS

1X LED FACE MASK

1X CONTROLLER

1X USB-C CABLE

1X USER MANUAL

IMAGES ARE FOR ILLUSTRATION PURPOSES ONLY. MADE IN CHINA.

FOR MORE INFORMATION OR TO DOWNLOAD MORE DETAILED DATA,
PLEASE CONTACT THE SELLER.

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WEAR AND TEAR OF THE PRODUCT AND ITS COMPONENTS CAUSED BY THEIR COMMON USE (DISCHARGED BATTERIES, LOSS OF THE ORIGINAL COLOR OF THE PRODUCT, DECREASED CAPACITY OF ACCUMULATOR, WEAR AND TEAR OF MECHANICAL OR DISPLAYING COMPONENTS, WEAR AND TEAR OF SUPPLY CABLES, ETC.) SHALL NOT BE SUBJECT NEITHER TO DEFECT WARRANTY NOR QUALITY WARRANTY IN THE MOMENT OF TAKE-OVER.
FURTHERMORE, COMPLAINTS ABOUT PRODUCT PARAMETERS, NOT INCLUDED IN THE OPERATING INSTRUCTIONS OR IN ANY COMMERCIAL-TECHNICAL DOCUMENTATION OF SELLER, RELATED TO THE PRODUCT, SHALL NOT BE SUBJECT TO QUALITY WARRANTY IN THE MOMENT OF TAKE-OVER.

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