# TERAPEUTICKÁ LED MASKA NA RTY

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## PALSAR7®

### PACKAGING INCLUDES:

- 1x LED Lip Mask
- 2x Mouthpiece
- 1x USB Cable
- 1x User Manual

### CHARGING

Connect the USB cable included in the package to the device. The indicator light will flash slowly while charging. When the device is fully charged, the indicator light will turn green. If you do not use the device for a long time, charge it at least once every 3 months to protect the battery.


The images with descriptions are for illustrative purposes only.

## THERAPEUTIC LED LIP MASK

This device uses red-light therapy to penetrate deep into the skin of the lips, promoting blood circulation and natural collagen production. Regular use helps improve the health of the lips from the inside out, contributing to a fuller appearance, more pronounced color and smoother texture. The device also helps minimize wrinkles, promotes volume and can reduce the damage caused by cold sores. Ideal for at-home lip care without the need for invasive procedures.

The device is equipped with 32 LEDs, which are made up of four separate light sources with different wavelengths in the range of 600-850 nm. This combination of lights helps to naturally fill the lips, unify their shade and smooth fine lines and wrinkles around the mouth.

## PRODUCT DESCRIPTION
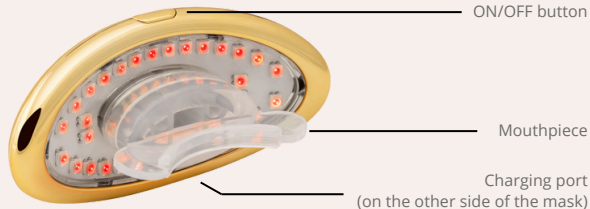
## RECOMMENDED USE

1. Wash and cleanse your lips thoroughly before use.
2. Attach the mouthpiece to the LED mask and gently bite down.
3. Hold down the button to turn on the device - a three-minute treatment cycle will automatically start.
4. Short press the button to switch between two levels of light intensity.
5. After 3 minutes, the device will automatically turn off.
6. Use a maximum of 3 times a day.

## LIGHT SENSITIVITY TEST

Before first use, perform a light sensitivity test.

1. Remove the mouthpiece from the device.
2. Turn on the LED mask and slowly move it over the skin on your forearm in circular motions for 3 minutes.
3. After application, observe the test area for at least 5 minutes.
4. If redness, irritation or any discomfort occurs, discontinue use.

## SPECIFICATIONS

- Rated voltage: 3.7V
- Rated power: 1.8W
- Charging voltage: DC-5V
- Charging current: 1A
- Product size: 86x60x58 mm
- Product weight: 61 g

## WARNINGS

- Before first use, perform a sensitivity test according to the procedure in this manual.
- Pay attention to cleanliness - always clean the device and the attachment thoroughly before and after use according to the instructions.
- Do not charge the device for longer than absolutely necessary.
- The product is not suitable for people with dental diseases, allergies, hypersensitivity to light (photosensitivity), cheilitis, pregnant women and children.
- If you suffer from a skin disease, consult a dermatologist before use.
- In case of any discomfort, stop using the device and contact a doctor.
- The device is not heat-resistant - protect it from fire and corrosive substances.
- Never immerse the mask in water and do not use water or chemicals to clean it - use only a slightly damp cloth.
- Do not expose the device to direct sunlight.
- Use only in accordance with the instructions for use and observe basic safety precautions.
- The product is for external use only.
- Do not bite or consume the product or any part of it.
- If the silicone tip is damaged, replace it with a spare one to avoid inhaling small parts.
- For hygiene reasons, it is not recommended to share the device with other people.
- Keep out of reach of children.
- Never use the device or charging cable if they are damaged.

## STORAGE AND MAINTENANCE

- After use, remove the mouthpiece from the device.
- Rinse the mouthpiece with water and let it air dry.
- Clean the device itself only with a damp cloth - never immerse it in water or rinse it under running water.
- Store the device in a dry and cool place, away from direct sunlight and moisture.

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SkyTech Corporation s.r.o., U libeňského pivovaru 2015/10,