GALVANIC EYE CARE MASSAGER

USER MANUAL

**Dear customer,
thank you for purchasing this
product.**

Please read this user manual carefully and keep it for future reference.

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EN 1-10

This manual contains the product's user manual in 4 languages.

For the manual in other languages using the QR code below.

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GALVANIC EYE CARE MASSAGER

PRODUCT INFORMATION

The Rechargeable Eye Massager based on the principle of magnetic electricity and ion infusion utilizes high-frequency micro vibration (13,000 times/min) and constant heating (20/40/42°C) technology to massage and lift skin around eyes, improve the absorption of skin care products, help reduce under-eye puffiness/eye bags/dark circles/crow's feet/fine lines/wrinkles. Meanwhile, the device is equipped with three kinds of beauty lights (blue, yellow and red) to effectively soothe skin, brighten skin tone and increase skin elasticity.

Consistent use of the device will bring you younger and refreshed-looking eyes.

WARNING

- Do not place or store it over or near water contained in a bathtub, washbasin, sink, etc.
- If you have a skin condition, please consult a dermatologist before use.
- Test the device on your arm prior to the first use.
- If you experience any discomfort, discontinue use immediately and consult a physician.
- The device is not waterproof. Do not rinse it under faucet or immerse it into water.

- For reasons of hygiene, we do not recommend sharing the device with anyone else.
- Keep the device out of reach of children to avoid accidents.
- Never use the device or charging cable if they are damaged.
- For external use only.
- Do not bite or ingest the product or any of its parts.
- Use only as directed and follow basic safety precautions.

PRODUCT OVERVIEW

PRODUCT FEATURES

- Constant Heating: 20/40/42°C temperature modes, meet different needs
- Micro Vibration: 13,000 times/min high frequency, lift and massage skin
- Beauty Light: blue, yellow and red beauty lights, care about your skin
- Fast Heating: reach the setting value in 15 seconds, wait for nothing
- Smart Memory: automatically select your last choice, more convenient
- Smart Induction: contact between hand and induction strip activates device
- Auto Shutdown: switch itself off after been left for 5 minutes

CHARGING METHODS

- The battery may run out during long distant transportation. It's recommended to charge the device prior to the first use.
- The charging indicator will turn red when the battery is running low. Please recharge in time.

1. Connect the device with an adapter (not included) via the Type-c charging cable.
2. Plug into the power supply. The charging indicator will keep red during charging.
3. The charging indicator will turn off automatically once the battery been fully charged. It takes about 2 hours.

MODE INTRODUCTION

- 20°C (Blue Light): Soothe and calm skin, eliminate eye bags
- 40°C (Yellow Light): Brighten skin tone, lighten dark circles
- 42°C (Red Light): Firm skin, fade crow's feet/fine lines/wrinkles

SPECIFICATIONS

Body Material: ABS

Battery Capacity: 600mAh

Charging Voltage: 5V

Charging Time: 2 Hours

Input Current: 270mA

Output Power: Max 2.5W

Product Size: 144.3*24.5*23.5mm

Product Weight: 52g

PACKAGE CONTENTS


1x PALSAR7 Rechargeable Eye Massager

1x Type-C Charging Cable

1x User Manual

EYE CONTOUR


1.

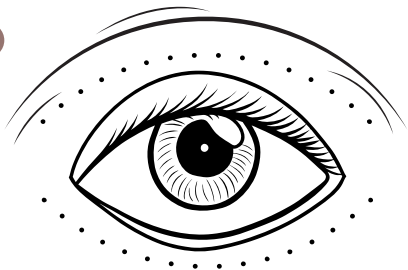Apply eye cream or eye serum to your eye contour.

Please apply proper amount of eye cream or serum. It's wasteful to apply too much eye cream or serum as they cannot be fully absorbed by the skin. What's more, the residual liquid may get inside of the device and shorten its service life.


2.

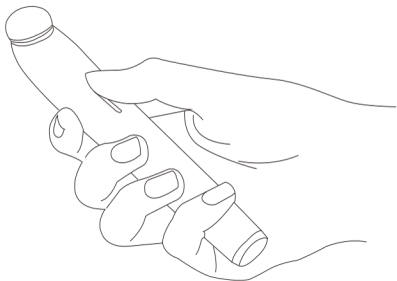Press and hold the Power/Mode Button for 2 seconds to turn the device on and short press to select a mode.

Note: The device will automatically select you last choice.


3.

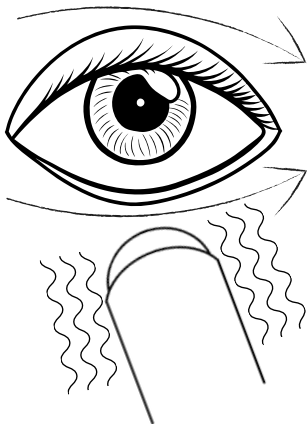Gently glide the massage head to lift and massage the eye contour. Repeat 15-20 times.

Make sure your hand directly contacts with the induction strip, otherwise the device will not work.

4.

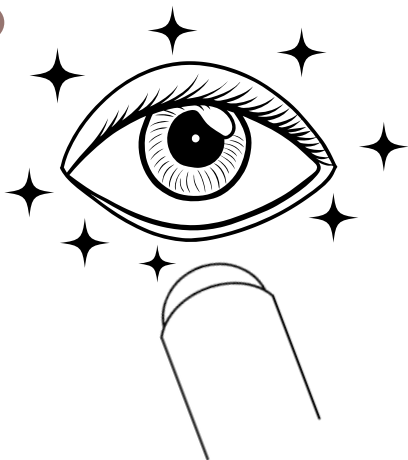Gently press the massage head to improve the absorption of eye cream or serum. Stay at each point for about 10 seconds and repeat 2-3 times.

5.

After use, press and hold the Power/Mode Button for 2 seconds to turn the device off. Then wipe the massage head with clean cloth.

WHOLE FACE


1.

Apply a facial care product (serum, essence, gel or cream).

Again, use an appropriate amount to avoid shortening the life of the device.


2.

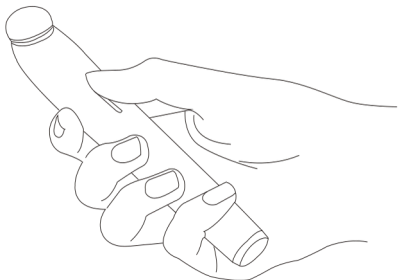Press and hold the Power/Mode Button for 2 seconds to turn the device on and choose a mode.

Note: The device will automatically select you last choice.

3.

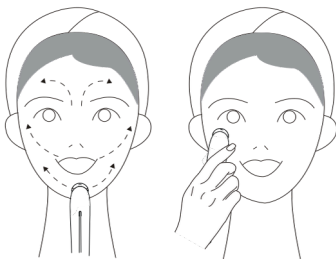Gently glide the massage head to lift and massage the whole face. Repeat 15-20 times. Help reduce forehead wrinkles/laugh lines/fine lines around mouth and accelerate blood circulation.

Make sure your hand directly contacts with the induction strip, otherwise the device will not work.

4.

After use, press and hold the Power/Mode Button for 2 seconds to turn the device off. Then wipe the massage head with clean cloth.

Note: The device can be used to massage other parts of your body, such as arms, waist and legs.

STORE & CARE

- Use clean cloth (dry or wet) to clean the device after each use.
- Do not use cleansers containing alcohol to clean the device.
- To protect the battery, recharge the device every 3 months when it is not in use for a long term.
- However, do not charge the device for longer than absolutely necessary.
- Do not leave the device in direct sunlight and never expose it to extreme heat.
- Store in a dry and shady place away from damp environments.

Images with descriptions are for illustration purposes only.

MADE IN CHINA.

For more information or to download more detailed data, contact your seller.

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Wear and tear of the product and its components caused by their common use (discharged batteries, loss of the original color of the product, decreased capacity of accumulator, wear and tear of mechanical or displaying components, wear and tear of supply cables, etc.) shall be subject neither to defect warranty nor quality warranty in the moment of take-over.

Furthermore, complaints about product parameters, not included in the operating instructions or in any commercial-technical documentation of Seller, related to the product, shall not be subject to quality warranty in the moment of take-over.

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