
PINK CLOUD

LED MASK 8IN1

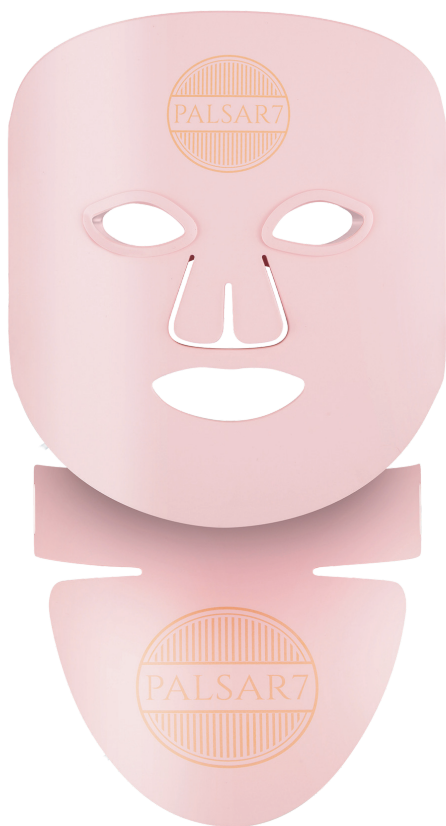PROFESSIONAL TREATMENT AT HOME

**DEAR CUSTOMER,
THANK YOU FOR PURCHASING THIS PRODUCT.**

PLEASE READ THIS USER MANUAL CAREFULLY AND KEEP IT FOR FUTURE REFERENCE.

THIS PRINTED MANUAL CONTAINS THE USER MANUAL IN 3 LANGUAGES. FOR OTHER LANGUAGE VERSIONS, PLEASE SCAN THE QR CODE BELOW.

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EN

DEAR CUSTOMER, THANK YOU FOR PURCHASING THIS PRODUCT. PLEASE READ THIS MANUAL CAREFULLY BEFORE USE AND RETAIN IT FOR FUTURE REFERENCE.

PALSAR7 - ABOUT THE BRAND

PALSAR7 IS A BRAND SPECIALIZING IN INNOVATIVE PRODUCTS FOR PROFESSIONAL SKIN CARE FROM THE COMFORT OF YOUR HOME. IT OFFERS MODERN TECHNOLOGIES THAT SUPPORT SKIN HEALTH AND BEAUTY WITHOUT THE NEED TO VISIT A BEAUTY SALON. THE BRAND EMPHASIZES QUALITY, SAFETY AND EASE OF USE OF ITS PRODUCTS. THANKS TO THIS, YOU CAN INDULGE IN PROFESSIONAL CARE ANYTIME, ANYWHERE.

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
EN ABOUT THE PRODUCT

THE DEVICE USES THE PRINCIPLE OF LIGHT THERAPY TO ACTIVATE SKIN CELLS AND IMPROVE SKIN CONDITION THROUGH SPECIALLY ADJUSTED LIGHT OF SPECIFIC WAVELENGTHS. IT IS A NON-INVASIVE, PAINLESS AND SAFE METHOD OF SKIN CARE WITHOUT THE USE OF CHEMICALS. AFTER BEING ABSORBED BY THE SKIN, THE LIGHT IS CONVERTED INTO INTRACELLULAR ENERGY, WHICH TRIGGERS A PHOTOCHEMICAL AND ENZYMATIC REACTION. DURING THIS PROCESS, A LARGE AMOUNT OF GLYCOGEN AND PROTEIN IS PRODUCED, WHICH PROMOTE CELL METABOLISM AND COLLAGEN REGENERATION. THANKS TO THIS, EFFECTS SUCH AS REJUVENATION, ACNE TREATMENT, DARK SPOT REDUCTION, SKIN FIRING, WRINKLE SMOOTHING AND SKIN WHITENING, CAN BE ACHIEVED.


PRODUCT DESCRIPTION

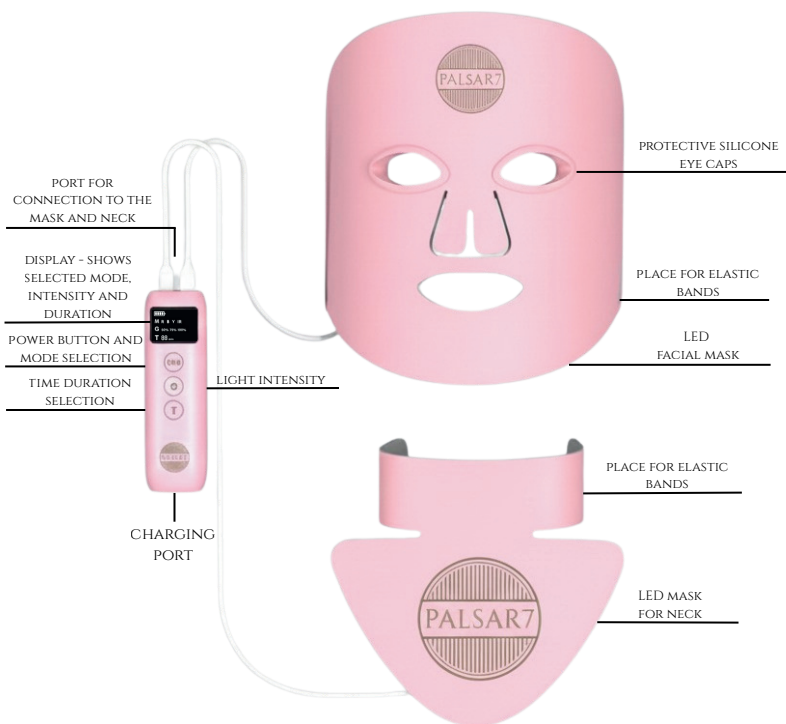
PRODUCT FUNCTIONS

THIS FACE LED MASK USES 72 PIECES AND THE NECK MASK USES 60 PIECES OF LEDS IN 4 COLORS (RED, INFRARED, GREEN, BLUE), THANKS TO THE COMBINATION OF LIGHTS, THE MASK OFFERS 8 MODES. THEIR FUNCTIONS ARE AS FOLLOWS:


RED LIGHT (630 NM \pm 10 NM): THIS LIGHT, CALLED BIOACTIVE LIGHT, IS CHARACTERIZED BY HIGH PURITY, INTENSITY AND UNIFORM DENSITY, WHICH MAKES IT SUITABLE FOR SKIN CARE AND HEALTH. RED LIGHT INCREASES CELL VITALITY, ACCELERATES BLOOD CIRCULATION, PROMOTES CELL METABOLISM, STIMULATES THE REGENERATION OF COLLAGEN AND FIBROUS TISSUE, THEREBY IMPROVING SKIN ELASTICITY AND REDUCING DULLNESS. IT IS A GREAT ANTI-AGING AID, ACTS AS AN ANTIOXIDANT AND PROMOTES SKIN RENEWAL.

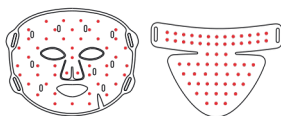
INFRARED LIGHT/INVISIBLE LIGHT (850 NM \pm 10 NM): THIS LIGHT PROMOTES COLLAGEN PRODUCTION IN THE SKIN, SPEEDS UP ITS METABOLISM, BRIGHTENS SKIN TONE, TIGHTENS PORES, AND HELPS ACHIEVE A MORE YOUTHFUL APPEARANCE. SINCE INFRARED LIGHT IS INVISIBLE, THE EMITTED LIGHT CANNOT BE SEEN WHEN YOU SELECT THIS MODE.

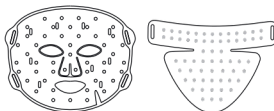
BLUE LIGHT (450 NM \pm 10 NM): THIS LIGHT CORRESPONDS TO THE ABSORPTION MAXIMUM OF PORPHYRIN, A METABOLITE OF PROPIONIBACTERIUM ACNES. DURING THE DE-EXCITATION PROCESS AFTER STIMULATION OF PORPHYRIN, A LARGE AMOUNT OF SINGLET ACTIVE OXYGEN IS GENERATED, WHICH CREATES A HIGHLY OXIDATIVE ENVIRONMENT THAT DESTROYS PROPIONIBACTERIUM ACNES. THIS HELPS IN THE TREATMENT OF ACNE.

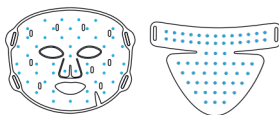
GREEN LIGHT (520NM \pm 10NM): IT CAN HELP SOOTHE THE SKIN, RELAX THE NERVES, AND RESTORE BALANCE TO SENSITIVE SKIN. STABILIZES AND RELAXES THE NERVOUS SYSTEM AND PROMOTES THE REMOVAL OF TOXINS FROM THE BODY.

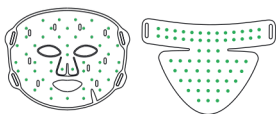
YELLOW LIGHT (COMBINATION OF RED AND GREEN LIGHT) (580 NM \pm 10 NM): THIS LIGHT HELPS TO IMPROVE THE OXYGEN EXCHANGE FUNCTION OF SKIN CELLS, PROMOTES MICROCIRCULATION, BREAKS DOWN AND LIGHTENS PIGMENT, STIMULATES CELL GROWTH, AND THUS SMOOTHS THE SKIN, BRIGHTENS SKIN TONE, REGENERATES THE SKIN, AND REDUCES WRINKLES.

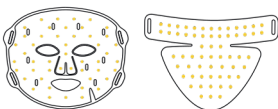
PURPLE LIGHT (COMBINATION OF RED AND BLUE LIGHTS) (420 NM \pm 10 NM): THIS LIGHT IS A GOOD HELPER IN TREATING ACNE AND REMOVING ACNE SCARS.

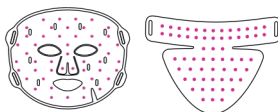CYAN LIGHT (COMBINATION OF BLUE AND GREEN LIGHTS) (490NM \pm 10NM): THIS LIGHT CAN INCREASE CELLULAR ENERGY AND PROMOTE METABOLISM.

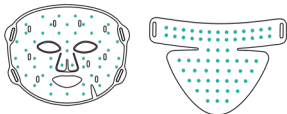WHITE LIGHT (COMBINATION OF RED, BLUE AND GREEN LIGHTS) (400-700 NM \pm 10 NM): THIS LIGHT CAN PENETRATE DEEP INTO THE SKIN, ACCELERATE METABOLISM, DECOMPOSE PIGMENTATION, AND IMPROVE FINE WRINKLES AND SKIN LAXITY.

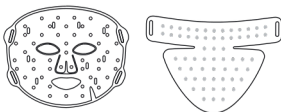HOW TO USE

1. WASH YOUR FACE AND NECK AND GENTLY DRY WITH A SOFT TOWEL. THEN APPLY TONER TO YOUR FACE AND NECK.


2. THEN YOU CAN APPLY YOUR FAVORITE NATURAL COSMETICS (SERUM, MASK OR ESSENCE).


3. TAKE OUT THE FACE AND NECK MASK. CONNECT THE USB-C CONNECTOR FROM THE FACE AND NECK MASK TO **THE TOP PORT** ON THE CONTROLLER.


4. HOLD THE MASK CLOSE TO YOUR FACE AND NECK. GRASP BOTH ENDS OF THE ELASTIC STRAPS AND FASTEN THEM WITH THE VELCRO AT THE BACK OF YOUR HEAD. ADJUST THE STRAPS TO ENSURE A SNUG FIT AND PREVENT LIGHT FROM ENTERING AROUND THE SILICONE EYE PROTECTION.

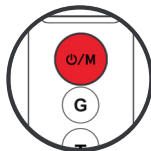5. PRESS AND HOLD THE POWER BUTTON FOR 1.5 SECONDS TO TURN ON THE DEVICE. AFTER POWERING ON, THE DEFAULT MODE IS RED LIGHT. SHORT PRESS THE MODE BUTTON TO SWITCH BETWEEN MODES. THERE ARE EIGHT MODES AVAILABLE: 1 - RED LIGHT (R), 2 - INFRARED LIGHT (IR), 3 - BLUE LIGHT (B), 4 - GREEN LIGHT (G), 5 - YELLOW LIGHT (Y), 6 - PURPLE LIGHT (R+B), 7 - CYAN LIGHT (B+G), 8 - WHITE LIGHT (R+B+G). YOU CAN ALSO SWITCH BETWEEN MODES DURING USE BY PRESSING THE MODE BUTTON.

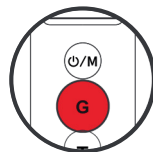6. SHORT PRESS TO ADJUST THE LIGHT INTENSITY. THERE ARE THREE BRIGHTNESS LEVELS: 50%, 75%, 100%.

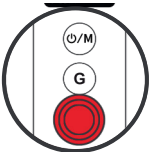7. SHORT PRESS THE TIME SETTING BUTTON TO SELECT THE PROCEDURE TIME. (5/10/15/20/25/30 MIN.)

8. AFTER YOU ARE DONE USING THE MASK, APPLY CREAM TO YOUR FACE AND NECK TO PREVENT DRYNESS.

RECOMMENDATION:
THE RECOMMENDED USAGE TIME IS 10-30 MINUTES 3-4 TIMES A WEEK.
USE THE BLUE LIGHT FOR A MAXIMUM OF 10 MINUTES PER TREATMENT.
DO NOT EXCEED THE MAXIMUM WEEKLY USAGE TIME (2 HOURS).

PRODUCT FEATURES

EFFECTIVE

HIGH-ENERGY, LARGE-AREA PHOTON SKINCARE WITH 4-IN-1 LEDS. FACE MASK: 72 LEDS; NECK MASK: 60 LEDS.

VERSATILE

8 LIGHT BEAUTY MODES TO ADDRESS VARIOUS SKIN CONCERNS.

SAFE

NON-INVASIVE SKIN CARE METHOD ENSURES CHEMICAL-FREE AND PAINLESS TREATMENT.

EASY TO USE

ONE-STEP THERAPY AND A BUILT-IN TIMER MAKE THIS MASK EASY TO USE.

COMFORTABLE

SOFT AND COMFORTABLE SILICONE MATERIAL ENSURES PLEASANT USER EXPERIENCE.

FLEXIBLE

THE MASK IS ONLY ATTACHED TO THE CONTROLLER, ALLOWING YOU TO MOVE FREELY DURING USE.

USABLE WITH COSMETICS

USE WITH COSMETIC PRODUCTS SUCH AS ESSENCE, SERUM OR MASK TO IMPROVE EFFECTIVENESS.

TIMING

6 TIMING MODES (5 MIN., 10 MIN., 15 MIN., 20 MIN., 25 MIN. AND 30 MIN.).

INTENSITY SETTINGS

3 ADJUSTABLE LIGHT INTENSITY LEVELS (50%, 75% AND 100%).

CLEAN

WIPEABLE, SILICONE DOES NOT TRANSFER BACTERIA, IS HYPOALLERGENIC AND EASY TO MAINTAIN.

CHARGING

- **PLEASE FULLY CHARGE THE DEVICE BEFORE USING IT FOR THE FIRST TIME.**
- CONNECT THE CONTROLLER TO THE APPROPRIATE ADAPTER (NOT INCLUDED) USING THE USB-C CABLE THAT WAS INCLUDED IN THE PACKAGE! **USE ONLY THE BOTTOM PORT FOR CHARGING!**
- PLUG INTO THE POWER OUTLET. WHILE CHARGING, THE SCREEN WILL LIGHT UP AND A BATTERY SYMBOL WILL APPEAR WITH THE CORRESPONDING CHARGE LEVEL.
- A FULL BATTERY SYMBOL INDICATES A FULL CHARGE. ONCE CHARGED, UNPLUG THE ADAPTER FROM THE POWER OUTLET AND YOU CAN START USING THE CONTROLLER AFTER CONNECTING IT TO THE MASK.

WARNING

- THE PRODUCT IS NOT SUITABLE FOR CHILDREN.
- DO NOT USE THE DEVICE IF YOU KNOW THAT YOU ARE PARTICULARLY SENSITIVE TO LIGHT AND HAVE NOT RECENTLY UNDERGONE A LIGHT SENSITIVITY TEST.
- DO NOT USE THE DEVICE IF YOU ARE PREGNANT AS THE RISKS ARE UNKNOWN AND HAVE NOT BEEN ESTABLISHED, CONSULT YOUR DOCTOR.
- DO NOT USE THE DEVICE IF YOU SUFFER FROM EPILEPSY AS SOME LIGHTS OR FREQUENCIES MAY TRIGGER A SEIZURE.
- THE DEVICE IS NOT WATERPROOF. DO NOT ALLOW WATER OR OTHER LIQUIDS TO ENTER ANY PART OF THE DEVICE. THIS MAY DAMAGE IT.
- DO NOT USE THE DEVICE IF YOU ARE TAKING MEDICATION THAT MAY MAKE YOUR SKIN SENSITIVE TO LIGHT.
- DO NOT USE THE DEVICE WHILE SLEEPING, WALKING OR DRIVING.
- USE THE DEVICE ACCORDING TO THE CORRECT INSTRUCTIONS AND AVOID EXCESSIVE USE.
- DISCONTINUE USE IF YOU EXPERIENCE DISCOMFORT OR IF REDNESS OR SKIN DISCOLORATION PERSISTS FOR MORE THAN 24 HOURS. CONSULT A DERMATOLOGIST.
- DO NOT TOUCH THE LEDs WHEN THE DEVICE IS TURNED ON.
- DO NOT USE THE DEVICE IF IT IS VISIBLY DAMAGED.
- DO NOT SUBJECT THE DEVICE TO SHOCKS TO AVOID DEFORMATION OR DAMAGE.
- DO NOT USE CHEMICALS, WATER, OR OTHER CLEANING OR LIQUID PRODUCTS FOR CLEANING.
- DO NOT EXCEED THE RECOMMENDED USAGE TIME.

HOW TO STORE

- AFTER USE, WIPE THE DEVICE WITH A DRY CLOTH. DO NOT RINSE UNDER WATER. DO NOT SOAK.
- STORE THE DEVICE IN A COOL, SHADY PLACE AND AVOID DIRECT SUNLIGHT OR HIGH TEMPERATURES.
- KEEP THE DEVICE OUT OF THE REACH OF CHILDREN TO PREVENT DAMAGE OR INJURY.
- NEVER ATTEMPT TO OPEN OR REPAIR THE DEVICE YOURSELF. IN THE EVENT OF A MALFUNCTION, CONTACT THE CUSTOMER SERVICE CENTER OR THE SELLER.
- AFTER USE, STORE THE DEVICE IN THE BOX.
- DO NOT PLACE ANY WEIGHT ON THE DEVICE. DO NOT PLACE IT UNDER OTHER OBJECTS AND BE CAREFUL TO AVOID MECHANICAL DAMAGE.
- WE RECOMMEND STORING IT IN A HORIZONTAL POSITION, INCLUDING THE CABLE. NEVER TRY TO DISCONNECT, DEFORM OR TWIST THE CABLE ATTACHED TO THE MASK.

SPECIFICATIONS

RATED CURRENT: C5V == 1A

RATED VOLTAGE: DC5V

RATED POWER: 6.5 W

CONTROLLER SIZE: 37X27X129 MM

CONTROLLER WEIGHT: 105 G

BATTERY CAPACITY: 5200MAH

FACE MASK SIZE: 284X217X20 MM

NECK MASK SIZE: 270X192X7 MM

MASK WEIGHT: 532 G

THE PACKAGE CONTAINS

1X LED FACE MASK

1X LED NECK MASK

1X ELASTIC BANDS FOR THE FACE AND NECK

2X SILICONE EYE PROTECTORS

1X CONTROLLER

1X USB-C CABLE

1X USER MANUAL

1X TEXTILE BAG

IMAGES ARE FOR ILLUSTRATION PURPOSES ONLY. MADE IN CHINA.

FOR MORE INFORMATION OR TO DOWNLOAD MORE DETAILED DATA, PLEASE CONTACT THE SELLER.

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WEAR AND TEAR OF THE PRODUCT AND ITS COMPONENTS CAUSED BY THEIR COMMON USE (DISCHARGED BATTERIES, LOSS OF THE ORIGINAL COLOR OF THE PRODUCT, DECREASED CAPACITY OF ACCUMULATOR, WEAR AND TEAR OF MECHANICAL OR DISPLAYING COMPONENTS, WEAR AND TEAR OF SUPPLY CABLES, ETC.) SHALL BE SUBJECT NEITHER TO DEFECT WARRANTY NOR QUALITY WARRANTY IN THE MOMENT OF TAKE-OVER. FURTHERMORE, COMPLAINTS ABOUT PRODUCT PARAMETERS, NOT INCLUDED IN THE OPERATING INSTRUCTIONS OR IN ANY COMMERCIAL-TECHNICAL DOCUMENTATION OF SELLER, RELATED TO THE PRODUCT, SHALL NOT BE SUBJECT TO QUALITY WARRANTY IN THE MOMENT OF TAKE-OVER.

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