# AquaFloss Station F300

---

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
## PRODUCT DESCRIPTION

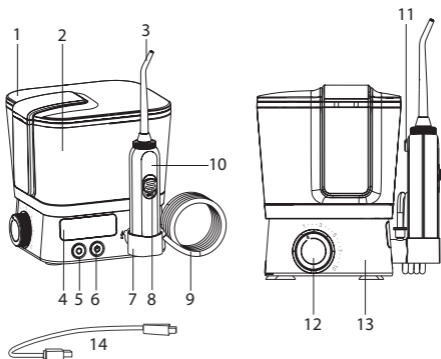1. Water tank lid/storage box
2. Water tank
3. Interchangeable head
4. Display
5. "Mode change" button
6. On/off button
7. Bracket
8. Switch
9. Connecting tube
10. Oral irrigator handle
11. Head release button
12. Water pressure setting rotary knob
13. Base
14. Charging cable

**Warning:** There may be residual water in the device. This is distilled water used during quality control testing.

## INSTRUCTIONS FOR USE

Before the first use, remove all packaging material.

If the water flosser has been stored in a low temperature environment during transport, leave it at room temperature for at least two hours before use.

1. Connect the head to the shower. A click indicates correct insertion. You can choose from the following headers:

- **Standard Head** – For daily plaque removal from interdental spaces and along gum lines
- **Dental Plaque Head** – For more effective removal of plaque around implants, crowns and bridges
- **Periodontal Head** – For cleaning periodontal pockets using a rubber tip
- **Orthodontic Head** – For removing plaque from brackets and wires of braces and other hard-to-reach places
- **Tongue Cleaner Head** – For effective cleaning of dirt on the tongue

2. Remove the water tank from the base. Remove the cap from the reservoir and fill the reservoir with water whose temperature does not exceed 40°C. The tank capacity is 600 ml. Then put the cap on. Make sure the water does not overflow.

3. Press the mode change button and select the desired shower mode. Modes:

- **Normal mode** – suitable for daily cleaning. Pressure Range: 20/35/68/109/125 PSI
- **Point mode** – suitable for effectively removing food residues between the teeth, thanks to a program of three quick and precise sprays under high pressure followed by a short pause. Pressure Range: 10/17/30/41/52 PSI
- **Pulse mode** – suitable for massaging teeth and gums. Pressure Range: 17/30/62/80/107 PSI

4. Use the rotary knob on the side of the base to set the pressure level (1-10 levels).

5. You can change the mode or pressure even while showering.

6. Leave your mouth slightly open to allow water to rinse from your mouth into the sink. Turn the switch on the handle to the ON position and press the On / Off button to start the flosser. Direct the water jet at a 90° angle to the teeth and gums. Slowly run the stream of water around your teeth. Direct the stream of water onto the gums and into the space between the teeth. If you have braces, you can direct the water between the joints of the braces.

7. When you have finished flossing, switch off the device using the On / Off button. Open the water tank cap and pour out any remaining water. Switch the device on again with the On / Off button to drain the water from the pump. Press the On / Off button to turn off the device and close the water tank cap.

**Warning:** To prevent the build-up of bacteria, ensure that no water remains in the device. If you do not intend to use the device for a longer time, make sure that it is thoroughly dried before storing it.

**Warning:** The oral shower must not be used for a long time with the outlet closed.

## CHARGING

The water flosser is not supplied with a fully charged battery. Before using for the first time, fully charge the water flosser. When charging for the first time, we recommend charging the oral shower for approximately 12 hours.

Remove the silicone cover of the charging port. Plug one end of the USB cable into the charging socket at the rear of the water flosser. Plug the other end of the USB cable into an adapter (adapter not included). Insert the adapter into the plug socket. The indicator light flashes red during charging. When the oral irrigator is charged, it will light up green. It takes approximately 6.5 to 8 hours to fully charge the water flosser.

**Warning:** To ensure maximum battery life, we recommend charging the water flosser once a week.

**Warning:** Charge the oral irrigator in a horizontal position on a stable surface. Keep the charging port dry and clean.

**Warning:** If you are not going to use the oral shower for a long time, we still recommend charging it regularly at least once every three months.

**Warning:** Do not use the oral shower while it is charging.

## INSERTING AND REPLACE THE NOZZLE HEAD

Insert the nozzle head into the hole at the top of the water flosser. You will hear a click when the nozzle head has been correctly attached.

To remove the nozzle head from the main body, press and hold the release button and then detach the nozzle head from the water flosser.

The nozzle head rotates 360°.


**Warning:** Only insert and remove the nozzle head when the water flosser is switched off.

**Warning:** Nozzle heads are items of personal hygiene and are not covered by warranty.

**Tip:** You can store the heads in a storage box.

## INSTALLATION OF THE HOLDER

Insert the end of the bracket into the hole in the front of the base at an angle of 15°. Press firmly and turn counter-clockwise until the bracket is straight.

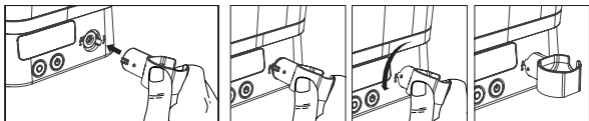Place the mouthwash in the holder after use.

To remove the holder, follow the opposite procedure.

## MAINTENANCE

Use only water and neutral detergents to clean the device. Never use caustic or abrasive cleaners as they may damage the device. Do not use water whose temperature exceeds 50 °C. Keep the device away from high temperatures and avoid direct sunlight.

**Warning:** Do not submerge the device in water.

### 1. The device

Remove dirt by wiping with a towel or cloth with liquid soap and water, then ensure all

water is rinsed off and wipe dry.  
Do not immerse the device in water.

**Warning:** Water must not enter the USB socket.

## 2. Water tank

Remove the water tank before cleaning. Wash the tank with water and wipe with a dry cloth.

Dry the tank thoroughly if you do not intend to use it for a longer time.

**Warning:** For greater ease of cleaning, the water tank can be removed from the body of the water flosser.

## 3. Nozzle head

Wash the nozzle head with water and wipe with a cloth.

## POSSIBLE COMPLICATIONS

Problem	Cause	Solution
The head cracks	The head has been used for a long time	Replace the head with a new one
The water tank is leaking	The seal on the body of the device is damaged	Contact the manufacturer with a request to provide a new seal
Water is leaking from the connection point of the head	The head is not attached correctly	Remove the head and reattach it when you hear a click
Water pressure is very low	Oral shower is used upside down	When using, the head must be placed higher than the water tank
	The water tank is empty	Fill the tank with water
	The water pressure level is set to low	Increase the water pressure level

## SPECIFICATIONS

Battery life	20 days (2x daily 2 minutes)
Modes	3 (Normal, Point, Pulse)
Pressure levels	10 (for each mode)
Timer	Yes (2 minutes)

Memory	Yes
Additional heads	F-series (Standard, Dental Plaque, Orthodontic, Periodontal, Tongue Cleaner)
Tank volume	600 ml
Charging	USB cable
Pressure range	30-130 PSI
Energy consumption	12W
Voltage	3.7V
Battery capacity	2000mAh
Water resistance	IPX5
Weight	570 g
Dimensions	130×135×140 mm

## PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

This product should be stored out of the reach of children. Use by children should take place only under the direct supervision of an adult.

## SAFETY GUIDELINES WHEN OPERATING THE DEVICE

1. Read the enclosed safety instructions carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: [www.truelife.eu](http://www.truelife.eu).
3. Printing errors and changes of information reserved.

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