InterdentBrush I3

PACKAGE CONTENTS

- 1× TrueLife InterdentBrush I3
- 30× Replacement heads
- 1× AA battery
- 1× Manual
- 1× Safety instructions

PRODUCT DESCRIPTION

- 1) Nylon cleaning head
- 2) Interdental brush head cover
- 3) Switch
 - 0: off
 - 1: low speed
 - 2: high speed
- 4) Handle
- 5) Battery cover


INSERTING THE BATTERY

The battery must be inserted into the interdental brush before use.

1. Carefully remove the battery cover and insert one battery. When inserting, observe the correct polarity, which is marked on the bottom of the compartment.
2. Put the cover back on.

Follow the same procedure when replacing the battery.

When the efficiency of the interdental brush decreases, the battery needs to be replaced.


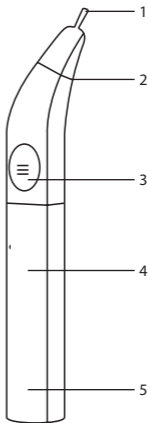
ATTACHMENT OF THE NYLON CLEANING HEAD

1. Rotate the toothbrush head cover until the dot symbol on the head is aligned with the  symbol on the handle.
2. Pull the cleaning head out of the handle. If it is attached too tightly, use the tool to remove the cleaning head on the storage box. Insert the cleaning

head all the way into the hole in the storage box and slide the handle together with the interdental brush head in the direction of the arrow. This releases and removes the cleaning head from the attachment.

3. Use the handle to pull the new cleaning head out of the storage box. Push the handle into the hole of one cleaning head and pull out the cleaning head after it has clicked into place.
4. Check that the cleaning head is firmly seated on the attachment. If necessary, push the cleaning head into the attachment with a little more force.
5. Replace the toothbrush head cover and rotate it so that the dot symbol on the head is aligned with the  symbol on the handle.

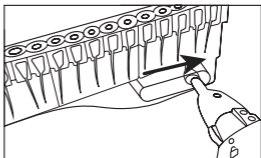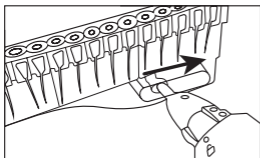CLEANING THE INTERDENTAL SPACE


We recommend using the interdental brush in front of a mirror to guide the cleaning head more confidently and prevent gum injuries.

The interdental toothbrush can be used every day to clean between the teeth. In this case, we recommend replacing the cleaning head approximately every 3 days, as the edges of the cleaning head may become bent, reducing the cleaning efficiency.

Keep the interdental brush in an upright position throughout the cleaning process.

We recommend dividing the teeth and interdental spaces into 4 sections. Once you have cleaned all the spaces between the teeth in one section, rinse the cleaning head under running water and move on to the next section.

1. Insert the cleaning head from the outside into the interdental space to be cleaned. The cleaning head is long enough so there is no need to clean the interdental spaces from behind.
2. Move the switch to position 1 for low speed or to position 2 for high speed. The cleaning head starts to move.
3. Move the interdental brush up and down to clean the entire lateral surfaces of adjacent teeth.

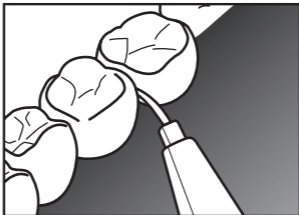Notice: Do not use force when using the interdental brush. If you experience discomfort during cleaning, select a low speed or move to another interdental space.

4. When the interdental space is cleaned, move the switch to position 0. The interdental toothbrush stops.
5. Pull the cleaning head out of the interdental space and continue cleaning through the next interdental space.

After cleaning, rinse the cleaning head thoroughly under running water. Hold the interdental brush so that the cleaning head is pointing downwards and the water flows over the cleaning head.

STORAGE

When not in use, store the toothbrush in a dry place away from children and pets. We recommend removing the battery from the toothbrush if you will not be using it for a long time.

SPECIFICATIONS OF THE INTERDENTAL TOOTHBRUSH

Battery life	75 days (2 times a day for 2 minutes), 5 hours
Number of movements	24,000 movements/min
Modes	2 (low speed and high speed)
Additional heads	I-series
Brush dimensions	190 × 23 × 55 mm
Weight of the brush	65 g
Waterproofing	IPX4
Power	1.5 V AA 1,000 mAh

MANUFACTURER RESPONSIBILITY

TrueLife products are warranted against manufacturing defects. In the event of such defects, elem6 s.r.o. is obliged to repair or replace the defective part or product with a comparable part or product at its discretion. Transfer, alteration or tampering with the product will exceed the scope of this warranty. This is the case even if the product is lost or damaged due to a defect, negligence or other error. The warranty does not cover damage resulting from use, accident or normal wear and tear. Elem6 s.r.o. assumes no liability for any accident, injury, death, loss or other claim related to or arising out of the use of this product. In no event shall elem6 s.r.o. be liable for incidental or consequential damages related to or resulting from the use of this product or any of its components. Returns or exchanges of replacement parts and/or products may be subject to shipping, handling, replacement and/or restocking charges.

This product should be stored out of the reach of children. Use by children is only allowed under adult supervision.

SAFETY INSTRUCTIONS:

1. Read the package safety information carefully.
2. For proper handling of the product, please read the current version of the manual, which is available for download at: www.truelife.eu.
3. Printing errors and information changes are subject to change.

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