**TESLA FLOSSER
FS200**
USER MANUAL

EN

DE

CZ

SK


PL

HU

SL

PARAMETRY PRODUKTU

| | |
|----------------------|--|
| Product Name | Water Flosser |
| Product Size | 66x49.5x166mm(Excluding nozzle height) |
| Reservoir Capacity | 200ML |
| Water Pressure | About 60~130psi |
| Frequency | 1400-1800times/min |
| Nozzle Diameter | 0.06mm |
| Flossing Mode | Standard(Spot-sparly), Soft, Strong, Pulse |
| Rated Input | 5 V, 1 A |
| Rated Power | 5W |
| Battery Capacity | 1000mAh |
| Charging Port | Type-C Charging Cable |
| Waterproof | IPX7 |
| Time of Long-lasting | About 30 days(1 minute per day based on Strong mode) |
| Accessories | Device - 1pc Nozzle - 4pcs User Manual - 1pc Type-C Charging Cable - 1pc |
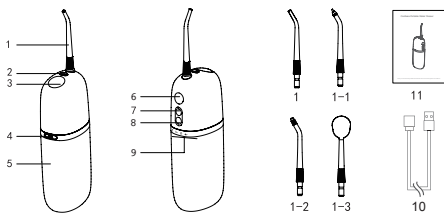
| Number | Description | Function |
|--------|-----------------------------|---|
| 1 | Standard Nozzle | Close to the gap between teeth, suitable for daily use. |
| 1-1 | Periodontal Pocket Nozzle | Close to the gap of the teeth, suitable for people with sensitive gums. |
| 1-2 | Orthodontic Nozzle | Can close to the gap between teeth to clean and washable braces, suitable for orthodontic people. |
| 1-3 | Tongue Cleaner Nozzle | Effectively remove tongue moss to avoid bacteria from entering the body. |
| 2 | Nozzle Release Button | Easy to load and unload the nozzle after pressing the button. |
| 3 | In-Reservoir Nozzle Storage | Insert the nozzle back into the handle. |
| 4 | Type-C Charging Port | Connect the charging device. |
| 5 | Reservoir | Store water. |
| 6 | Power/Flossing Button | Start the water spray / Off the water spray. |
| 7 | Pressure Control Button | Adjust the water spraying gear, blue low-grade, white mid-grade, orange high-grade. |
| 8 | Mode Button | Adjust the point of spray, soft, strong, pulse. |

| | | |
|----|-----------------------|---|
| 9 | Mode Indicator | Select the mode corresponding indicator lights up. |
| 10 | Type-C Charging Cable | Connect the product and power adapter for charging. |
| 11 | User Manual | Product instructions of guiding users. |

CHARGING

Charging The Water Flosser

Plug the Type-C charging into the charging port of water flosser, and USB connector to a power adapter with 5V output(No included).

Suggest to charge over 3 hours if meet below situation:

- (1)It is the first time to use this device.
- (2)The device is to be stored for an extended period of time (3 months or more).
- (3)The device is low in power or stops working during use.

Note:

- Please put the metal objects such as coins or clips away from USB charging cable while charging, as this may cause the device to overheat.
- It is normal with heating when charging. If still shorten working after fully charged, it maybe at the end of the battery life. Please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.
- Use this device at the temperature between 0 °C~35 °C.

USING METHOD

Kindly note that turn Water Flosser off and charge over 3 hours before first use.

1. Insert nozzle into the center of the knob at the top of the water flosser handle. You will hear a „click” once the nozzle is correctly locked into. (Suggest to replace the nozzle every 3 months for better results. If the nozzle is loose or blocked by foreign objects, it should be replaced).
2. Disconnect charging cable from the unit prior to filling the reservoir and product use. Lift flip top on reservoir or take the water tank, and fill the reservoir with lukewarm water.
3. Close the nozzle collection cover or rotate and locked the water tank to prevent the water tank from overflowing.
4. Shortly press the mode button “⊞”, to switch the flossing mode.

Note: The device will keep same mode with last shutdown for next working due to the memory function.

- Standard Spray Mode: Holding the power button to spray water, flossing after 1 second, and release the button to stop spraying water.
- ☞ Soft Mode: For the users with sensitive teeth or first time users.
- ⊞ Strong Mode: For the normal cleaning of teeth with most of peoples.
- ⊞ Pulse Mode: Alternating strong and weak, deep cleaning, and massaging gums.
- ≡ Pressure Control: Water pressure adjustment. Low in Blue light, Medium in White light, Strong in Orange light.

5. Shortly press the pressure control button “≡”, to switch the intensity of the water spray.
6. Ensure the nozzle kept parallel to your teeth and gums. Direct the jet stream at a 90-degree angle to your gumline. Slightly close lips to avoid splashing but allow the water to flow out into your sink.
7. Shortly press the power/flossing button, start flossing with






white lighting. Support to switch gears and modes when working. The longest flossing time is 90 seconds then turn off automatically.

8. Aim the tip toward teeth, glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

METHOD OF USE

Note:

The device can be turned on with Power Button/Flossing Mode Button/Pressure Control Button. Select the required mode and water pressure for cleaning. And then start flossing after pressing the Power/Flossing Button. The device will automatically turn off after 5 seconds if no operation is performed.

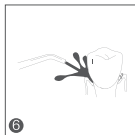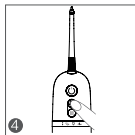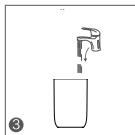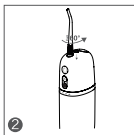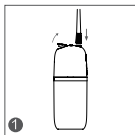
CLEANING & CARE

Do not use any corrosive or abrasive detergent as it will damage the water flosser.

Cleaning

1. Press and hold the nozzle release button to remove the nozzle, rinsed with water. Please clean every time after used.
2. The water tank can be detached, clean it with a soft cloth to wipe exterior and inside.
 - Please do the clean as above steps and store in a dry environment if you plan not to use this product for a long time.

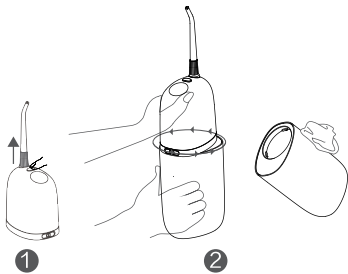Battery Care and Disposal

1. At the end of the battery life, please recycle the battery at your nearest recycling center per local or state requirements. If you can not remove the battery, please take the water flosser to the specified recycle center.
2. Do not remove the battery by users. It may be damaged in the process of forcibly removing the battery, so please operate with caution.

TROUBLESHOOTING

| Problem | Possible Cause | Solution |
|--|--|--|
| Do not work. | (a) New one, bought just now. (b) Stored over 3 months. | (a) Charge the water flosser over 3 hours. (b) Contact service. |
| The water flosser only operates for several minutes after charged. | (a) No fully charged, low-battery. (b) At the end of battery lifetime. | (a) Charge the water flosser over 3 hours. (b) Contact service. |
| The water pressure is low. | (a) The nozzle is blocked. (b) The nozzle is deformed. (c) The pump is blocked. | (a) Replace the nozzle. (b) Contact service. |
| No flossing happened. | (a) There is no water in the reservoir. (b) There is not correct angle when using the flosser. (c) The water pipe is not assembled well. (d) The pump is blocked. | (a) Fill water into the reservoir. (b) Direct the nozzle at 90-degree angle to gumline. (c) Assembled the water pipe well. (d) Contact service. |

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed.

DANGER:

To reduce the risk of electrocution:

- Always unplug the product after using.
- Do not handle the plug with wet hands.
- Do not place in or drop into water or any other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a bathtub, or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.

- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:

Make sure you turn the Water Flosser off and charge for more than 3 hours before the first use.

Note:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device.
- If the product is opened/disassembled for ANY reason, the warranty is VOID.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if dropped into any water.

Contact the manufacturer or your local distributor.

- Do not directly wash under the tongue, in the ear, nose or other delicate areas. This product is capable of producing pressure that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by a dental professional.
- Only use tips and accessories recommended by the manufacturer.
- Do not drop or insert any foreign object into any opening or hose.
- Keep the electrical cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

IMPORTANT SAFEGUARDS

- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. They can reduce product performance and will shorten the life of the product. Use of these products may void the warranty.
- Remove any oral jewellery prior to use of this product.

- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Do not use this product for more than 5 minutes in each two-hour period.

NOTE:

Closely instruct and supervise children 8 years and above (children younger than 8 years should not use the product) and individuals with special needs in the proper use of this product.

- Children should be supervised to ensure that they do not play with the appliance.
- Fill the water tank only with warm water or other solution recommended by a dental professional.

INFORMATION ABOUT DISPOSAL AND RECYCLING

This product is marked with the symbol for separate collection. The product must be disposed of in accordance with the regulations for disposal of electrical and electronic equipment (Directive 2012/19/EU on waste electrical and electronic equipment). Disposal together with regular municipal waste is prohibited. Dispose of all electrical and electronic products in accordance with all local and European regulations at the designated collection points which hold the appropriate authorisation and certification in line with the local and legislative regulations. Correct disposal and recycling help to minimise impacts on the environment and human health. Further information regarding disposal can be obtained from the vendor, authorised service centre or local authorities.

EU DECLARATION OF CONFORMITY

Hereby, Tesla Global Limited declares that the radio equipment type TSL-PC-FS200 is in compliance with EU directives. The full text of the EU declaration of conformity is available at the following internet address: tsl.sh/doc

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