# ToothWave

H7001

- DE BEDIENUNGSANLEITUNG
- EN USER MANUAL
- FR MODE D'EMPLOI
- BG РЪКОВОДСТВО ЗА УПОТРЕБА
- HU HASZNÁLATI UTASÍTÁS
- RO INSTRUCȚIUNI DE UTILIZARE
- CS NÁVOD K POUŽITÍ
- SL PRIROČNIK ZA UPORABNIKA
- PL INSTRUKCJA OBSŁUGI

Patent protected  
technology by

# TOOTH WAVE™


SCIENCE DRIVEN ORAL CARE

DentalRF™  
Technology

# Silk'n®

# ToothWave

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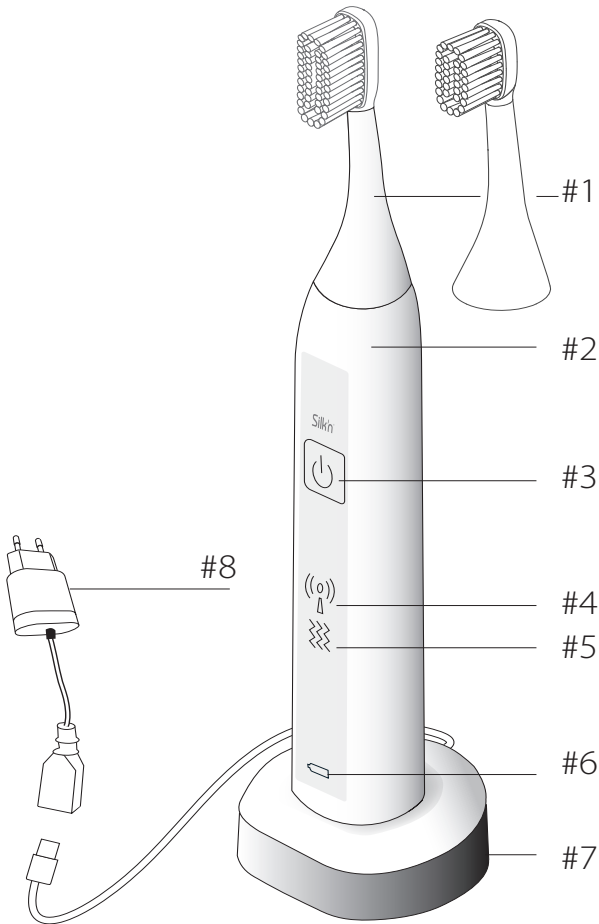DE

# Table of contents

|  |           |
|--|-----------|
| <b>1. Warnings &amp; contraindications.....</b>      | <b>21</b> |
| <b>2. Signs &amp; symbols .....</b>                  | <b>25</b> |
| <b>3. Get to know Silk'n ToothWave .....</b>         | <b>26</b> |
| 3.1 Intended use.....                                | 26        |
| 3.2 Package contents.....                            | 26        |
| 3.3 Device controls & settings.....                  | 27        |
| 3.4 Silk'n ToothWave's brush head .....              | 28        |
| 3.5 Device charging.....                             | 28        |
| 3.6 Reading & storing the user manual.....           | 29        |
| <b>4. Treatment plan.....</b>                        | <b>29</b> |
| <b>5. What to expect with Silk'n ToothWave .....</b> | <b>29</b> |
| <b>6. Get started!.....</b>                          | <b>30</b> |
| <b>7. Cleaning &amp; maintenance .....</b>           | <b>30</b> |
| 7.1 Cleaning the toothbrush.....                     | 30        |
| 7.2 Cleaning the charging cradle.....                | 31        |
| 7.3 Replacing the brush head .....                   | 32        |
| <b>8. Parts &amp; consumables .....</b>              | <b>32</b> |
| <b>9. Troubleshooting .....</b>                      | <b>32</b> |
| <b>10. Customer service .....</b>                    | <b>33</b> |
| <b>11. Warranty .....</b>                            | <b>33</b> |
| <b>12. Disposal .....</b>                            | <b>33</b> |
| 12.1 Battery removal .....                           | 33        |
| 12.2 Packaging.....                                  | 35        |
| 12.3 Device .....                                    | 35        |
| <b>13. Specifications .....</b>                      | <b>36</b> |
| <b>14. Legal notice .....</b>                        | <b>37</b> |

## 1. Warnings & contraindications

---

A statement or event that indicates a possible or impending danger, problem, or other unpleasant situation when using the device.

### **Electricity & safety warnings**

- ⦿ Always ensure that your hands are dry before connecting the charging cradle to the electricity supply.
- ⦿ Do not place or store the USB wall adapter where it can fall, or be pulled into any liquid.
- ⦿ For the purpose of recharging the battery, only use the USB wall adapter provided with this device.
- ⦿ Do not use the USB wall adapter if it has fallen into liquid.
- ⦿ If the charging cradle malfunctions, unplug immediately.
- ⦿ Never use the charging cradle or the USB wall adapter if they are damaged or aren't working properly.

- ⦿ Keep the device, all parts and accessories away from heated surfaces.
- ⦿ Do not attempt to modify, open or repair the device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void your warranty.
- ⦿ Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if they are damaged in any way, do not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device, and contact customer service.
- ⦿ The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- ⦿ Do not clean any of the parts in the dishwasher.
- ⦿ Children shall not play with the device. Cleaning and user maintenance shall not be made by children.
- ⦿ Children must not play with the packaging. They may injure themselves with the packaging material, or swallow small parts and choke.

- ⦿ This device contains batteries that are non-replaceable.
- ⦿ Do not dispose of the device or the battery with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
- ⦿ If the product reaches its end of life, the battery must be removed and disposed of safely, before disposing of the product (see section 12. Disposal).

## **Contraindications**

This device including all its parts must not be used by the following people: persons below the age of 18, adults with limited physical, sensory or psychological capacities, those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.

## **Do not use if you:**

- have a pacemaker, internal defibrillator, another active implanted device, or if you have any medical concerns about using this device.
- have a current or past history of oral cancer or oropharyngeal cancer, or if you have any other type of cancer, or pre-malignant lesions.
- have an impaired immune system due to immuno-suppressive diseases, or you use immuno-suppressive medication.

## **Consult your physician before use if you:**

- have concurrent conditions, such as epilepsy, or cardiac disorders.
- are pregnant or nursing.

## **Consult your dentist before use if you:**

- have had oral or gum surgery in the last two months.
- have severely injured gums.

## **Additional warnings:**

- Do not bite the brush head.

## 2. Signs & symbols

|   |   |
|---|---|
|  | Important information about using or maintaining the device.  |
|  | Products marked with this symbol fulfill the requirements of the EU directives.   |
| IP67  | Protected from immersion in water with a depth of up to 1 meter (or 3.3 feet) for up to 30 mins.  |
|  | Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.                                  |
| <b>Li-ion</b>   | Lithium-ion battery inside.   |
|  | The built-in rechargeable ToothWave battery contains substances that may pollute the environment. Always remove the battery before you discard at an official collection point for batteries. |
|  | Products marked with this symbol fulfill the Korean safety requirements for electrical and electronic equipment.  |

## 3. Get to know Silk'n ToothWave

---

### 3.1 Intended use

Silk'n ToothWave is intended to reduce tartar (calculus\*), plaque\*, tooth discoloration\*, gingivitis\* and gum bleeding\*. Good oral care can prevent the development of tooth decay. The device is not a substitute for regular visits to the dentist or hygienist for routine dental care. See our websites and videos for more information about our patented DentalRF™ technology.

Silk'n ToothWave is for personal home-use only, and is not intended for use on multiple patients in a dental practice or institution. Do not use for any purpose other than as described in this manual.

\* After 6 weeks, ToothWave showed a statistically significant reduction of tartar (calculus), stains, gum inflammation and gum bleeding compared to the control toothbrush - SmileSonic Pro Advanced Clean.

### 3.2 Package contents

Refer to the diagram on page (i) at the front of the manual.

- 2x Brush heads **(#1)**
- 1x Handle **(#2)**
- 1x Charging cradle with USB cable + transparent base **(#7)**
- 1x USB wall adapter **(#8)**
- 1x Warranty booklet **(#9)**
- 1x User manual **(#10)**

### 3.3 Device controls & settings


Refer to the diagram on page (i) at the front of the manual.

| Control  | Function / Mode   |
|--|---|
| Control button<br><b>(#3)</b><br>                 | <b>Long press:</b> turns the device on or off<br><b>Short press:</b> each press moves to the next mode  |
| DentalRF mode indicator light<br><b>(#4)</b><br>  | DentalRF is active for all modes<br>   |
| Vibration mode indicator light<br><b>(#5)</b><br> | <b>Mode 1:</b> Low vibration   |
|  | <b>Mode 2:</b> Medium vibration    |
|  | <b>Mode 3:</b> High vibration<br>Recommended mode                                      |
|  | <b>Mode 4:</b> No vibration    |
| Battery indicator light<br><b>(#6)</b><br>      | <b>Green blinks</b> - battery is charging<br><b>No light</b> - battery is charged<br><b>Orange blinks</b> - battery is low<br><b>Orange constant</b> - battery is empty |

- i** When the battery is fully charged, the light will be constant green for 30 seconds, then it will go out.

### 3.4 Silk'n ToothWave's brush head

The brush head (#1) is an innovative and special feature of Silk'n ToothWave. It contains bristles, 2 DentalRF electrodes and a silicone strip that lies between the electrodes. The silicone strip enables the DentalRF current to flow across the teeth. It also ensures that the toothpaste and the DentalRF energy interact effectively. For this reason, we recommend that the bristles are always in good contact with the teeth.

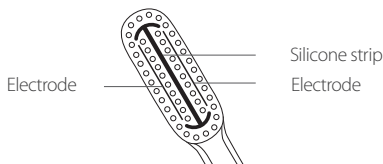
⚠ Do not use the brush head (#1) if the bristles, silicone strip, or electrodes become loose, bent, damaged, or crushed. Damaged electrodes, bristles or a damaged silicone strip may break off while brushing.

### 3.5 Device charging

- Check that the battery is not low before each use.

First charge   
hours

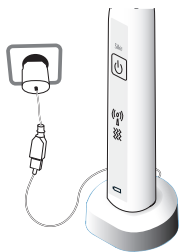
### 3.6 Reading & storing the user manual

This user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Read the manual in full before use, and refer to it if you are in any doubt about how to use the device safely.

## 4. Treatment plan

- Use for 2 minutes, morning and evening.
- Clean between the teeth and gently massage gums.
- Brush each section (1,2,3,4) for 30 seconds.
- i** The device is suitable for use with dental implants.
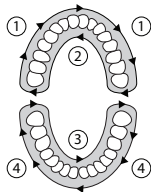
## 5. What to expect with Silk'n ToothWave

During or after using Silk'n ToothWave you may:

- feel a slight sensation of warmth. Don't worry, this is normal.
- feel pain, if you have cavities or any teeth condition that requires dental treatment. If you do feel pain, consult your dentist.
- i** Stop using the device if you experience excessive bleeding after use, or if bleeding continues after one week of regular use.

## 6. Get started!

---

1. Rinse the bristles under lukewarm water.
2. For best results, apply a generous amount of toothpaste. You can use any commercial toothpaste.
3. Turn on the device. Short press the control button **(#3)** as many times as needed to select your treatment mode.
4. Move the brush in circular movements, slowly and systematically from tooth to tooth, and gently massage the gums. Make sure the bristles are always in good contact with your teeth. Silk'n ToothWave comes with soft bristles to protect your teeth and gums. Brush gently, and apply only light pressure when brushing. Do not scrub.
5. Brush for 2 minutes. You will feel a slight vibration every 30 seconds (1 vibration at 30 seconds, 2 vibrations at 60 seconds and so on). After the vibration(s), you can move to another section of the mouth (refer to section [4. Treatment plan](#)). After 2 minutes, the device will shut down automatically.
6. Clean the device and store in open air whilst it is wet.

## 7. Cleaning & maintenance

---



### 7.1 Cleaning the toothbrush

It is important to clean the device after each use. Use lukewarm water.

- Clean the brush head **(#1)** and handle **(#2)**, especially if your toothpaste contains baking soda, or other bicarbonates. This is to prevent the plastic from cracking.


- Once a month, remove the brush head (#1) from the handle (#2). Clean the groove in the handle (#2) to remove any toothpaste or debris that may have accumulated. Use a cotton bud, rinse and dry with a clean dry cloth.
- Clean the base of the brush head (#1) in the same way.

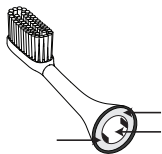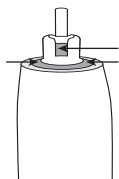EN

## 7.2 Cleaning the charging cradle


1. First disconnect the charging cradle (#7) from the electricity supply and remove the toothbrush from the charging cradle (#7).
2. Next, remove the transparent base. To do this, turn the charging cradle (#7) upside down. Remove the transparent base by pressing the clips outwards.

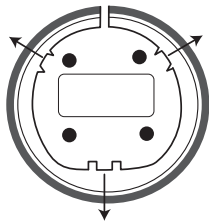3. Clean both parts with a damp cloth and dry them. Place the charging cradle (#7) back, and press down firmly. You will hear it click into place.

## 7.3 Replacing the brush head

Replace the brush head (**#1**) every 3 months, or whenever you see that the bristles are damaged or worn. To change the brush head (**#1**), simply pull out the brush head (**#1**) and replace with a new one.
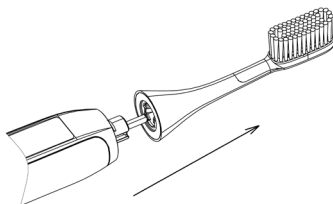
## 8. Parts & consumables

---

You can purchase replacements from where you purchased the device, or online at [www.silkn.eu](http://www.silkn.eu).

## 9. Troubleshooting

---

| Problem   | Checks  |
|---|---|
| My device does not start.                         | Make sure the battery is not low. Refer to section <a href="#">3.5 Device charging</a> .  |
| The device is charged but it does not work.       | Press the control button ( <b>#3</b> ) to restart the device. If the problem persists, contact customer service.                          |
| All indicator lights ( <b>#4, #5, #6</b> ) blink. | There is a system error. Press the control button ( <b>#3</b> ) to restart the device. If the problem persists, contact customer service. |

**i** For problems or queries, please contact Silk'n customer service.

## 10. Customer service

For more information see our website: [www.silkn.eu](http://www.silkn.eu). This manual is also available as a PDF download from the website. Contact Silk'n customer service to report unexpected operations or events, or any other problem with your device. Our contact details are:

| Country         | Service number     | Email                          |
|-----------------|--------------------|--------------------------------|
| Belgium         | 0800-29316         | serviceconsommateurbe@silkn.eu |
| France          | 0891-655557        | serviceconsommateurfr@silkn.eu |
| Germany         | 089 51 23 44 23    | kundenservicede@silkn.eu       |
| Italy           | +31 (0)180-330 550 | servizioclientiit@silkn.eu     |
| Netherlands     | 0900-2502217       | klantenservicenl@silkn.eu      |
| Spain           | 900 823 302        | servicioalcliente@silkn.eu     |
| UK              | 0906-2130009       | customerareuk@silkn.eu         |
| Other countries | +31 (0)180-330 550 | info@silkn.eu                  |

## 11. Warranty

Please refer to the separate product warranty booklet for complete information.

## 12. Disposal


### 12.1 Battery removal

At the product's end of life, you are required to remove the battery before disposal. To do this, the device must be opened up. The device will no longer function after this procedure, nor can it be repaired. Any warranty will be void when the device casing is opened and the battery is removed. Home Skinovations Ltd is not liable for any physical or material damages as a result of performing this procedure.

1. Remove the device from the charging cradle (#7) and make sure the battery is empty.


2. Remove the brush head (#1) from the handle (#2).

3. **Open the casing:** Use a small tool (e.g., screwdriver or utility knife), to detach the bottom cover from the rest of the device.

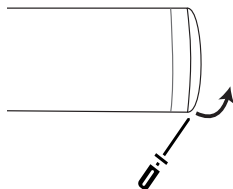
4. **Remove the inner section:**

Place the handle (#2) upside down on a sturdy table. Use the tool to move the clips inwards. As the clips release, you will be able to remove the inner section and lay it on the table.

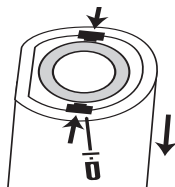5. **Remove the battery:** Use cutters or small scissors, to cut each of the 3 wires separately (red, black, yellow). Make sure the red (positive) and black (negative) wires do not touch each other. Then remove the battery carefully. Do not damage the battery.
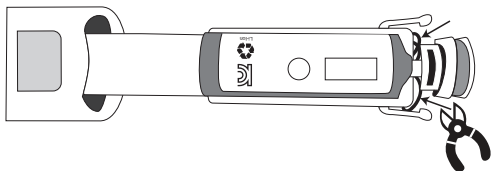
## 12.2 Packaging

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

## 12.3 Device

Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances\* into the environment. Remove the battery and dispose of it safely before disposing of the device.

\* Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

## 13. Specifications

|  |  |
|--|--|
| Model no.  | H7001  |
| Technology   | DentalRF™ and vibration  |
| Vibration  | 0; 275Hz; 300Hz or 400Hz (± 5%)  |
| Radio frequency  | 3.3MHz ±10%; 3W (max)  |
| Package size   | (W)165 (H)227 (D)80 [mm]   |
| System weight  | 115g   |
| Transport & storage between uses and storage condition | Temperature: - 40 to +70°C<br>Relative humidity: 10 to 90%rH<br>Atmospheric pressure: 500 to 1060hPa |
| Operating conditions                                   | Temperature: 5 to 40°C<br>Relative humidity: 15 to 90%rH<br>Atmospheric pressure: 700 to 1060hPa     |
| Mode of operation                                      | Non-continuous   |
| Adapter  | HX-S4V0500600 (Europe)   |
| Input  | 100 – 240V; 50/60Hz; 0.2A  |
| Output   | 5.0Vdc; 0.6A   |
| Average active efficiency                              | 73.12%   |
| No-load power consumption                              | 0.064W   |
| Commercial Registration No.                            | 91330281X10172885G   |
| Manufacturer   | Yuyao Hongxiang<br>Electrical Appliance Factory  |
| Address  | Xiaodong Industrial Zone, Yuyao,<br>Zhejiang Province, P.R.China                                     |

|                   |  |
|-------------------|--|
| Waterproof rating |  |
| Adapter           | IPX4 - Protects from splashing water, no matter the direction.                     |
| Charging cradle   | IPX7 - Protects in water up to 3 feet (1 meter).                                   |
| Handle            | IP67 - Totally protects against dust and protects in water up to 3 feet (1 meter). |

## 14. Legal notice

---

Copyright 2020 © Home Skinovations Ltd. All rights reserved.

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification. Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O.Box 533, Yokneam 2069206, ISRAEL.

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Copyright 2020 © Home Skinovations Ltd. Alle Rechte vorbehalten.
Copyright 2020 © Home Skinovations Ltd. Tous droits réservés.
Copyright 2020 © Home Skinovations Ltd. Toate drepturile rezervate.
Copyright 2020 © Home Skinovations Ltd. Všechna práva vyhrazena.
Copyright 2020 © Home Skinovations Ltd. Wszelkie prawa zastrzeżone.
Invention Works B.V.