# Silk'n®

Beautiful Technology

## SONICSMILE (DELUXE)

SG-986

EN USER MANUAL

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PT MANUAL DO UTILIZADOR

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DA BRUGERMANUAL

## **Legal notice**

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks,

copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd.

Invention Works B.V.


Donk 1B, 2991 LE Barendrecht

The Netherlands

[www.silkn.eu](http://www.silkn.eu) – [info@silkn.eu](mailto:info@silkn.eu)

# Table of contents

|                                                                     |    |
|---------------------------------------------------------------------|----|
| Overview.....                                                       | 5  |
| Package contents/device parts.....                                  | 6  |
| General information.....                                            | 7  |
| Reading and storing the operating instructions.....                 | 7  |
| Explanation of symbols.....                                         | 8  |
| Safety.....                                                         | 9  |
| Proper use.....                                                     | 9  |
| Safety instructions .....                                           | 10 |
| First use.....                                                      | 20 |
| Checking the toothbrush and package contents.....                   | 20 |
| Basic cleaning.....                                                 | 21 |
| Operation.....                                                      | 21 |
| Charging the rechargeable battery .....                             | 21 |
| Charging with the travel charging case (Deluxe version).....        | 22 |
| Using the toothbrush.....                                           | 23 |
| Brushing your teeth .....                                           | 24 |
| Brushing your tongue with the tongue cleaner (Deluxe version) ..... | 26 |
| Selecting the mode .....                                            | 26 |
| Cleaning.....                                                       | 28 |
| Storage.....                                                        | 29 |
| Troubleshooting .....                                               | 30 |
| Technical data.....                                                 | 32 |
| Declaration of conformity .....                                     | 34 |
| Disposal .....                                                      | 34 |
| Disposing of the packaging .....                                    | 34 |
| Disposing of the toothbrush .....                                   | 34 |
| Warranty.....                                                       | 36 |
| Customer service.....                                               | 37 |

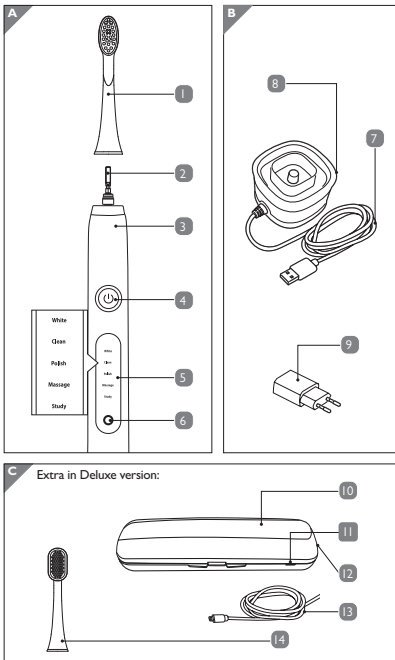\*Note: The content of the box may vary based on the model purchased.

## **Package contents/ device parts**

- 1 Brush head, 2×
- 2 Brush attachment point
- 3 Hand piece
- 4 Multifunction button
- 5 LED mode indicator
- 6 Charge level indicator
- 7 Main cord with USB port
- 8 Charging station
- 9 USB adapter
- 10 Travel charging case
- 11 Charge level indicator
- 12 Micro USB port
- 13 Micro USB cable
- 14 Tongue cleaner, 2x

## General information

### Reading and storing the operating instructions

These operating instructions accompany this SonicSmile (De-luxe) SG-986 (referred to below only as the “toothbrush”). It contains important information about start-up and handling. Before using the toothbrush, please read through the operating instructions carefully. This particularly applies to the safety instructions. Failure to follow these operating instructions may result in severe injury or damage to the toothbrush. The operating instructions are based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the operating instructions for future use. If you pass the toothbrush on to third parties, please be absolutely sure to include these operating instructions.

## **Explanation of symbols**

The following symbols and signal words are used in these operating instructions, on the toothbrush or on the packaging.

**WARNING!**

This signal symbol/word designates a hazard with moderate degree of risk which may lead to death or severe injury if not avoided.

**CAUTION!**

This signal symbol/word designates a hazard with low risk that, if not avoided, may result in minor or moderate injury.

**NOTICE!**

This signal word warns against potential damages to property.

This symbol provides you with useful supplementary information on assembly or operation.

Declaration of conformity (see chapter “Declaration of conformity”): Products marked with this symbol meet all applicable Community regulations for the European Economic Area.

## **Safety**

### **Proper use**

The toothbrush is designed exclusively for cleaning teeth. It is only intended for private use and is not suitable for commercial purposes.

Only use the toothbrush as described in these operating instructions. Any other use is deemed improper and may result

in damage to property or even personal injury. The toothbrush is not a children's toy.

Silk'n accepts no liability for damage caused by improper or incorrect use.

## **Safety instructions**

**WARNING!**

### **Risk of electric shock!**

A faulty electrical installation or excessive mains voltage may result in an electric shock.

- Only connect the charging station or charging case if the mains voltage of the socket corresponds to the specification on the nameplate.
- The main cord of the charging station cannot be replaced. In the

event of damage, the charging station must be disposed of.

- Only connect the charging station or charging case to an easily accessible socket so that you can quickly disconnect it from the mains in case of a failure.
- Do not operate the toothbrush, charging station or charging case if they exhibit visible damages or the main cord or USB port is defective.
- Do not open the housing; instead, have a qualified professional perform the repairs. Contact a qualified workshop for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection of the device or incorrect operation.

- Only parts that comply with the original device information may be used for repairs. Electrical and mechanical parts, which are essential for providing protection against sources of danger, are located in this toothbrush.
- The charging station or charging case must not be operated with an external timer or with a separate telecontrol system.
- Do not immerse the charging station, main cord, or charging case in water or other liquids.
- Never touch the USB port with damp hands.
- Never pull the USB port out of the socket by the main cord.
- Never use the main cord as a

carrying strap.

- Keep the charging station or charging case, toothbrush and USB port away from open flames and hot surfaces.
- Do not kink the main cord and do not lay it over sharp edges.
- Only use the charging station or charging case in rooms located indoors.
- Never set up the charging station or charging case such that there is a risk of it falling in a tub or sink.
- Never attempt to grab the charging station or charging case if it has fallen into water. In such a case, immediately disconnect the USB port.
- Only charge the toothbrush

with the enclosed charging station or charging case. Only use the USB adapter and micro USB cable provided.

- Make sure that children do not insert any objects in the charging station, charging case, or hand piece.
- Disconnect the USB port when you are not using the charging station or charging case, after the charging process, before you clean the toothbrush or put it away or in the event of a fault.

**WARNING!**

### **Risk of chemical burns!**

Fluid leaked from the rechargeable battery may lead to chemical burns.

- Do not allow acid from the rechargeable battery to come into contact with skin, eyes or mucous membranes.
- In the event of contact with acid from the rechargeable battery, immediately flush the applicable areas with plenty of clean water and consult a physician.

**WARNING!****Danger of explosion!**

Improperly handling rechargeable batteries may cause an explosion.

- Never short circuit the rechargeable batteries.
- Never dismantle rechargeable batteries.

- Never throw rechargeable batteries in fire.
- Never expose rechargeable batteries to excessive heat or direct sunlight.
- Do not in any case open the rechargeable batteries.

## **WARNING!**

**Danger for children and persons with impaired physical, sensory or mental capacities (e.g. partially disabled persons, older persons with reduced physical and mental capacities) or lack of experience and knowledge (e.g. older children).**

- This toothbrush may be used by children and persons with

impaired physical, sensory or mental capacities or those lacking experience and knowledge if they are supervised or have been instructed in how to safely use the toothbrush and have understood the risks associated with operating it. Children must not play with the toothbrush.

- Cleaning and user maintenance must not be performed by children.
- Keep children under eight years of age away from the toothbrush.
- Do not leave the toothbrush unattended while in use.
- Do not allow children to play with the small parts. Children could swallow small parts and choke.

## **CAUTION!**

### **Risk of injury!**

Handling the toothbrush improperly may result in injury.

- Lay the main cord so that it does not pose a tripping hazard.
- Immediately replace worn or damaged brush heads to prevent injury to the mouth and tooth area.

## **NOTICE!**

### **Risk of damage!**

Improper handling of the toothbrush may damage it.

- Do not drop the hand piece.  
If this does happen, do not

use the toothbrush any more. Have it checked by an authorised service point.

- Turn the toothbrush off if you would like to change the brush head.
- Place the charging station or charging case on an easily accessible, level and dry surface.
- Ensure that the main cord, the charging station or charging case and the toothbrush do not come into contact with hot parts.
- Stop using the toothbrush if its plastic parts exhibit cracks or breaks or are deformed.

## First use

### Checking the toothbrush and package contents

#### NOTICE!

#### Risk of damage!

If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the toothbrush.

- Be very careful when opening the packaging.
- 1. Take the toothbrush and its accessories out of the packaging.
- 2. Check to make sure that the delivery is complete (see **Fig. A, B** and **C** (the latter for the Deluxe version)).
- 3. Check whether the toothbrush or the

individual parts are damaged. If this is the case, do not use the toothbrush. Please contact the Silk'n service centre in your area (see chapter "Customer Service").

## **Basic cleaning**

1. Remove the packaging material and all protective foils.
2. Clean all parts of the toothbrush before first use as described in the chapter "Cleaning".

## **Operation**

### **Charging the rechargeable battery**

1. Insert the main cord with USB port **7** into the included USB adapter **9**.
2. Place the toothbrush on the charging station **8**. The toothbrush will now be inductively charged.

3. The charge level indicator **6** is illuminated in red during the charging process.
4. The toothbrush is fully charged after approx. 16 hours and the charge level indicator goes out.  
Once fully charged, the rechargeable battery will keep for up to 3 weeks.

## **Charging with the travel charging case (Deluxe version)**

1. Place the hand piece **3** in the travel charging case **10**. Insert the micro USB cable **13** into the micro USB port **12** of the case and the other end into the USB adapter **9**. Then, insert the USB adapter into a wall socket. The toothbrush will now be inductively charged.  
You can also charge via any CE certified device with a USB port. In this case, insert the micro USB cable into the micro USB

port of the case and the other end into the USB port of the device.

2. The charge level indicators **6** **||** on the toothbrush and case are illuminated in red during the charging process.
3. The toothbrush is fully charged after approx. 16 hours and the charge level indicators go out. Once fully charged, the rechargeable battery will keep for up to 3 weeks.

## Using the toothbrush

**CAUTION!**

### **Risk of injury!**

Handling the toothbrush improperly may result in injury.

- Do not apply too much pressure to the toothbrush and your teeth when brushing them.

- Your gums may bleed the first few times you use the toothbrush. If your gums bleed for more than 14 days, please consult a dentist.

- Brush your teeth as close to your gums as possible.
- Brush your teeth using small circular movements.
- Clean up to the spaces between your teeth.

## Brushing your teeth

1. Attach a brush head **1** to the hand piece **3** so that the brush attachment point **2** of the hand piece extends into the brush head.

The brush head must audibly and noticeably lock into place.

The brush head must point towards the multifunction button **4**.

2. Hold the brush head in clean water for a moment to wet it.

3. Apply a pea-sized amount of toothpaste to the brush head.
4. Carefully move the brush head up to your teeth.
5. Select the desired mode by pushing the multifunction button within 3 seconds (see chapter “Selecting the mode”).  
Wait for more than 3 seconds to switch the toothbrush off.
6. Brush the upper and lower parts of your teeth as well as the exterior and interior of your teeth for 30 seconds each without applying pressure.

After every 30 second interval, the toothbrush will stop for 0.5 seconds and automatically switch off after 2 minutes.

7. After each use, clean the toothbrush as described in the chapter “Cleaning”.

## Brushing your tongue with the tongue cleaner (Deluxe version)

1. Wet the tongue cleaner **14** and brush your tongue with the micro bristle side in a back and forth motion for 30 seconds. Clean as far back as comfortable.
2. Rinse your mouth and use the scraping side on your tongue for 30 seconds.

## Selecting the mode

The LED mode indicator **5** displays the current mode.

| Function | Operation | Description                                                |
|----------|-----------|------------------------------------------------------------|
| White    | Push 1 ×  | For white teeth and less staining.                         |
| Clean    | Push 2 ×  | For removing plaque, perfect for sensitive teeth and gums. |

| Function | Operation | Description                                                                                                                                   |
|----------|-----------|-----------------------------------------------------------------------------------------------------------------------------------------------|
| Polish   | Push 3 ×  | For polishing teeth.                                                                                                                          |
| Massage  | Push 4 ×  | For gently massaging your gums and improving the health of your gums.                                                                         |
| Study    | Push 5 ×  | Setting for getting used to the toothbrush. The vibration intensity gradually increases from weak to strong during the first 12 applications. |

## Cleaning

### NOTICE!

#### Risk of damage!

Handling the toothbrush improperly may result in damage.

- Do not in any case attempt to clean the toothbrush, charging station or charging case in the dishwasher. These parts would be destroyed as a result.
- Do not submerge the charging station or charging case in water or other liquids.
  1. After every use, hold the brush head  under running water and switch the toothbrush on.
  2. Remove all leftover toothpaste from the brush head and the hand piece  and

- switch the toothbrush back off.
3. Pull the brush head off of the hand piece.
  4. Also clean the brush head by holding it under running water.
  5. Dry the brush head and the hand piece.  
Clean the charging base and/or charging case with a soft, damp cloth and wipe it dry.
  6. Change the brush head no later than every 3 months.

As soon as the colour of the brush head fades, change it. You can purchase new brush heads on our website at [www.silkn.eu](http://www.silkn.eu).

## Storage

All parts must be completely dry before being stored.

- Always store the toothbrush in a dry area.
- Use the original packaging, a similar box or the travel charging case for storage.

- Store the toothbrush so that it is securely locked away preventing access by children.

## Troubleshooting

| Problem                        | Possible cause                                      | Solution                                                                                  |
|--------------------------------|-----------------------------------------------------|-------------------------------------------------------------------------------------------|
| The toothbrush isn't charging. | The charging station or charging case is defective. | Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours. |
|                                | The USB connector jack is defective.                | Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours. |
|                                | The charging station or charging case is defective. | Replace the charging station or charging case within the warranty period.                 |

| <b>Problem</b>                                                 | <b>Possible cause</b>                  | <b>Solution</b>                                                                                                                                   |
|----------------------------------------------------------------|----------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|
| The toothbrush operates with only a small amount of vibration. | The re-chargeable battery is depleted. | Charge the toothbrush for at least 16 hours.                                                                                                      |
| The toothbrush does not run for 2 minutes.                     | The re-chargeable battery is depleted. | Charge the toothbrush for at least 16 hours and make sure that the toothbrush is properly placed on the charging station or in the charging case. |

## Technical data

Model: SG-986

### Charging station

Input: DC 5V

Power consumption: 1.5 W

Protection class: II

IP code: IPX7

Dimensions

(L × W × H): 6.2 × 6.2 × 4.3 cm

### Hand piece

Rechargeable battery: 3.7 V, 700 mAh

Operating time: 42 times/full charge

Timer: 0.5 second pause  
after approx. 30  
seconds, deactivation  
after 2 minutes

Dimensions

(L × W × H): 2.6 × 2.7 × 25.4 cm

**Adapter**

|                   |                               |
|-------------------|-------------------------------|
| Model number:     | HNFG050060UE                  |
| Input:            | 100–240 V~,<br>50/60 Hz, 0.2A |
| Output:           | 5.0V $\equiv$ 600mA           |
| Protection class: | II                            |

**Travel charging case (Deluxe version)**

|                            |                      |
|----------------------------|----------------------|
| Model number:              | SG-209               |
| Working voltage:           | DC 5V $\equiv$ 300mA |
| Power rating:              | 1.5W                 |
| Dimensions<br>(L × W × H): | 22.6 × 6.4 × 3.6 cm  |

## Declaration of conformity

The EU declaration of conformity can be requested from the address specified on the back.

## Disposal

### Disposing of the packaging

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

### Disposing of the toothbrush

(Applicable in the European Union and other European states with systems for the separate collection of reusable waste materials)

## Old devices must not be disposed of with household waste!

If at some point it is no longer possible to use the toothbrush, each consumer **is required by law to dispose of old devices separately from their household waste**, e.g.

at a collection point in their community/borough. This ensures that old devices are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

## Batteries and rechargeable batteries do not belong in the household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances\* or not, to a collection point run by the communal authority or borough or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire toothbrush (with the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

\*labelled with: Cd = cadmium, Hg = mercury,  
Pb = lead

## **Warranty**

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## Customer service

For more information about the products of Silk'n visit your regional Silk'n website: [www.silkn.eu](http://www.silkn.eu). If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center.

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© 2020 Home Skinovations Ltd. All rights reserved. SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd.
InnoEssentials Deutschland GmbH