
PHILIPS


Lumea

BRI973, BRI976,
BRI977

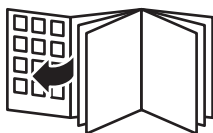1

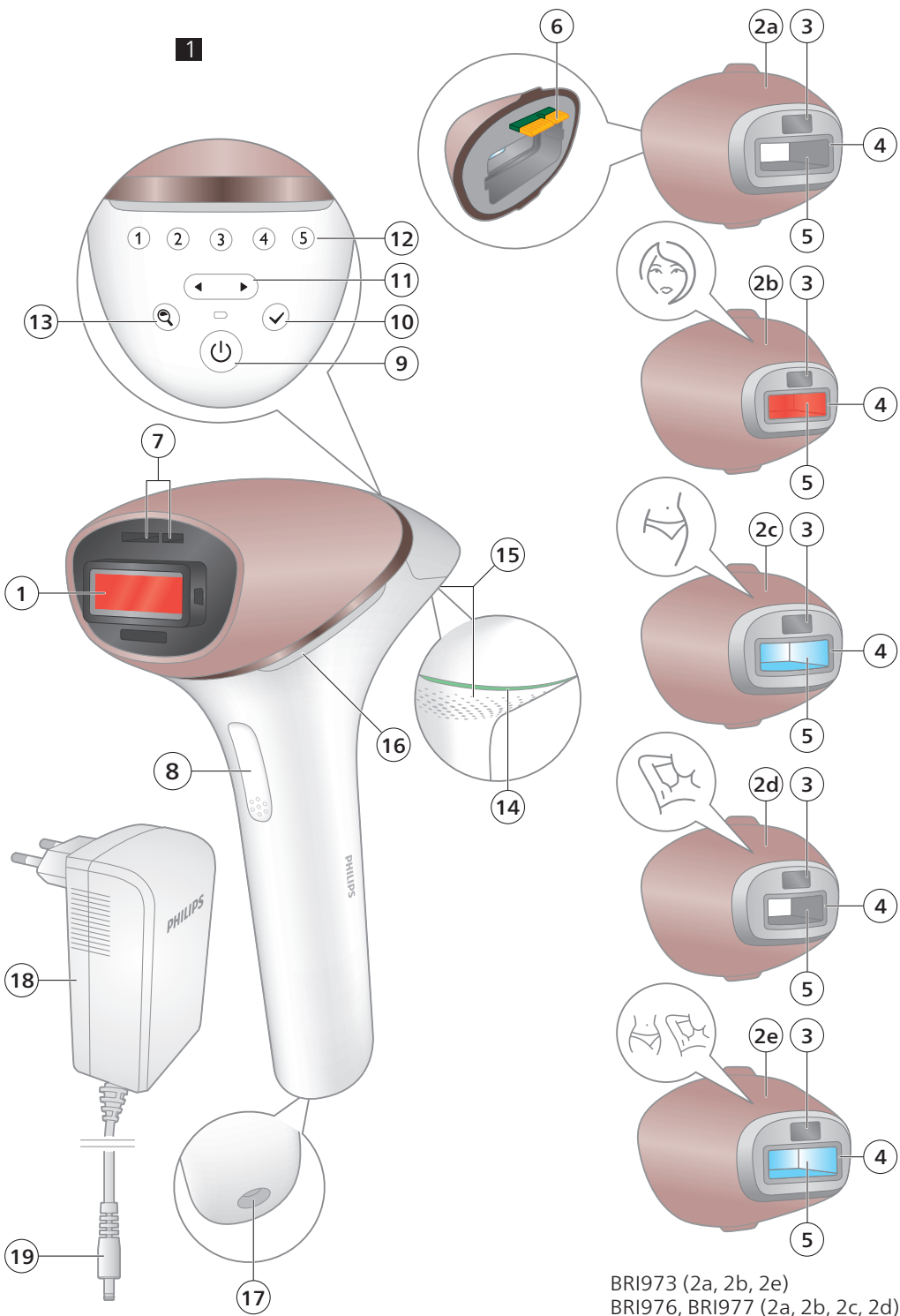BRI973 (2a, 2b, 2e)
BRI976, BRI977 (2a, 2b, 2c, 2d)

Suitable skin tones / Vhodné odstíny pleti /
Příkladne nijanse kože / Megfelelő bőrtónusok /
Odpowiednie odcienie skóry / Tonuri de piele
potrivite / Primerni odtenki kože / Vhodné
odtiene pokožky / Подходящи тонове на кожата

Suitable body hair colors / Vhodné barvy
chloupků / Pogodne boje dlačica / Alkalmas
testszőrzetszínek / Odpowiednie kolory włosów
na ciele / Culori adecvate ale părului de pe corp /
Primerne barve dlačic / Vhodné farby telesného
ochlpenia / Подходящ цвят на телесно
окосяване

English 6
Čeština 31
Hrvatski 57
Magyar 83
Polski 109
Română 137
Slovenščina 163
Slovensky 189
Български 215

Contents

Welcome	6
Device overview	7
Contraindications	7
Important safety instructions	9
How IPL works	11
Suitable body hair colors	12
Recommended treatment schedule	12
What to expect	13
How to use your Lumea Series 9900	14
Learning more about your Lumea Series 9900	19
Charging	19
Attachments	20
SmartSkin sensor	21
Selecting the light intensity	22
Two treatment modes: Stamp & Flash and Slide & Flash	22
Using your Lumea Series 9900 before and after tanning	23
Traveling with the device	23
After use	24
Troubleshooting	25
Warranty and support	27
Recycling	27
Technical specifications	28
Electromagnetic compatibility - Compliance information	29
Explanation of symbols	29

Welcome

Welcome to the beauty world of Philips Lumea Series 9900! You are only a few weeks away from silky-smooth skin.

Lumea Series 9900 uses Intense Pulsed Light (IPL) technology, known as one of the most effective methods to continuously prevent hair growth. In close collaboration with expert scientists and dermatologists, we adapted this light-based technology, originally used in professional beauty salons, for easy and effective use in the comfort of your home. Lumea Series 9900 is gentle and offers convenient and effective treatment at a light intensity that you find comfortable. Unwanted hairs are finally a thing of the past.

Lumea Series 9900 has been designed and developed for women, but may be suitable for use by men.

To fully benefit from the support that Philips offers, register your product at www.philips.com/support. For further information, please go to www.philips.com/Lumea to find our experts' advice, tutorial videos and FAQs and make the most of your Lumea.

Note: Keep these instructions with your product at all times.

Device overview

- | | |
|--|--------------------------------|
| 1 Light exit window with integrated UV filter | 8 Flash button |
| 2 Attachments | 9 On/off button |
| a Body attachment | 10 Confirmation button |
| b Facial attachment  | 11 Toggle buttons |
| c Bikini attachment  (BRI976, BRI977) | 12 Intensity lights |
| d Underarm attachment  (BRI976, BRI977) | 13 SmartSkin button |
| e Precision attachment  (BRI973) | 14 'READY to flash' light |
| 3 SmartSkin sensor | 15 Air input vents for the fan |
| 4 Integrated safety system | 16 Air output vents |
| 5 Reflector inside the attachment | 17 Device socket |
| 6 Electronic contacts | 18 Adapter |
| 7 Opening for electronic contacts | 19 Small plug |
| | 20 Luxurious pouch (not shown) |
| | 21 Cleaning cloth (not shown) |

Contraindications

Conditions

Never use the device:

- If you are under the age of 18 years or above 65 years.
- If you have skin type VI (you rarely to never sunburn, very dark tanning). In this case you run a high risk of developing skin reactions, such as hyperpigmentation and hypopigmentation, strong redness or burns.
- If you are pregnant or breastfeeding as the device was not tested on pregnant or breastfeeding women.
- If you have any active implants such as a pacemaker, neurostimulator, insulin pump, etc.

Medications/History

Never use the device if you take any of the medications listed below:

- If your skin is currently being treated with or has recently been treated in the past week with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical tretinoin and azelaic acid.
- If you have taken any form of isotretinoin Accutane or Roaccutane in the last six months. This treatment can make skin more susceptible to tears, wounds and irritations.
- If you are taking photosensitizing agents or medications, check the package insert of your medicine and never use the device if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.

- If you take anticoagulation medications, including heavy use of aspirin, in a manner which does not allow for a minimum 1-week washout period prior to each treatment.

Never use the device:

- If you have received radiation therapy or chemotherapy within the past 3 months.
- If you are on painkillers which reduce the skin's sensitivity to heat.
- If you take immunosuppressive medications.

Pathologies/Disorders

Never use the device:

- If you have diabetes or other systemic or metabolic diseases.
- If you have congestive heart disease.
- If you have a disease related to photosensitivity, such as polymorphic light eruption (PMLE), solar urticaria, porphyria, etc.
- If you have a history of collagen disorder, including a history of keloid scar formation or a history of poor wound healing.
- If you have epilepsy with flashlight sensitivity.
- If your skin is sensitive to light and easily develops a rash or an allergic reaction.
- If you have a skin disease such as active skin cancer, you have a history of skin cancer or any other localized cancer in the areas to be treated.
- If you have a history of vascular disorder, such as the presence of varicose veins or vascular ectasia in the areas to be treated.
- If you have any bleeding disorder.
- If you have a history of immunosuppressive disease (including HIV infection or AIDS).
- If you have a skin pigmentation disorder such as vitiligo.

Skin conditions

Never use the device:

- If you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex (cold sores), wounds or lesions and haematomas (bruises) in the areas to be treated.
- On irritated (red or cut), sunburned, recently tanned or fake-tanned skin.
- On the following areas: on moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. This can result in a burn and a change in skin color, which makes it potentially harder to identify skin-related diseases.
- On the following areas: warts, tattoos or permanent make-up.

Location/Areas

Never use the device on the following areas:

- Around the eyes and on or near the eyebrows. Use the device only below the cheekbones.
- On the lips, scalp, perineum, nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears.
- If you have had surgery in the areas to be treated.
- Men must not use the device on the face and neck including all beard-growing areas, nor on the whole genital area.
- On areas where you use long-lasting deodorants. This can result in skin reactions.
- Over or near anything artificial like silicone implants, subcutaneous injection ports (insulin dispenser) or piercings.

Note: This list is not exhaustive. If you are not sure whether you can use the device, we advise you to consult your doctor.

Important safety instructions

Warnings


Please read and follow these instructions to optimize the use of the device and prevent/minimize the risk of injury, skin reactions and side effects.

- Before you use Lumea, you should clean your skin and make sure it is hair-free, entirely dry and free from any lotion or gel.
- Keep the device and the adapter dry.
- The device is not washable. Never immerse the device in water and do not rinse it under the tap.
- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. near a filled bath, a running shower or a filled swimming pool). Moisture can damage the device and lead to a potential safety hazard.
- Always check the device before you use it. Do not use the device or adapter if it is damaged. Always replace a damaged part with one of the original type.
- If the device or adapter is broken or damaged, do not touch any inner part to avoid electric shock.
- Do not use the device if the UV filter of the light exit window and/or attachment is broken, as safe operation cannot be guaranteed under these circumstances.
- Do not attempt to open or repair your device or change the non-replaceable batteries. Opening Philips Lumea may expose you to dangerous electrical components and to pulsed light energy, either of which may cause bodily and/or eye injury.


- Always return the device to a service center authorized by Philips for examination or repair. Repair by unqualified people could cause an extremely hazardous situation for the user.
- No modification of this equipment is allowed. Do not modify, override or cut off any part of the device (e.g. battery, adapter or the cord), as this causes a hazardous situation.
- Do not look at the flash while using the device. Use the device in a well-lit room so that the light is less glaring to your eyes. While not necessary for safety, you may wear sunglasses or goggles for your convenience if you find the light uncomfortable or too bright.
- Do not flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase efficacy.
- For hygienic reasons, the device should only be used by one person.
- Do not use hair removal creams in the area (to be) treated with the device, as chemicals may cause skin reactions.
- Do not use any pencil or pen to mark the areas to be treated. This may cause burns on your skin.
- The device can be stored at a temperature between -25 °C and 70 °C. Charge and use the device at a temperature between 5 °C and 35 °C. The device will not work properly outside of this temperature range. This may cause burns on your skin.
- Only use this product for its intended purpose and follow the general and battery safety instructions as described in this user manual. Any misuse can cause electric shock, burns, fire and other hazards or injuries.
- Hair removal by intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.
- If you experience skin irritation or skin redness after treatment, wait until it disappears before applying any product to your skin. If you experience skin irritation after applying a product to your skin, wash it off with water. If you still get skin reactions, stop using the device and consult your doctor.
- If shaving, epilating or waxing causes skin irritation, we advise you not to use the device until the skin irritation has disappeared.
- This device is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the device. Supervision is required to avoid potential risks such as exposure to the light output, electrocution, strangulation from the cables or choking on the attachments.
- Tanning with natural or artificial sunlight might influence the sensitivity and color of your skin. Perform a skin test to determine the appropriate light intensity setting.
- After each treatment, wait at least 48 hours before tanning to ensure any skin reaction has been resolved.

- If you notice a skin tone change since the last treatment, we advise you to perform a skin test and to wait 24 hours before your next treatment.
- Always check the device before you use it. Make sure that the SmartSkin sensor window is clean when using the device, to avoid improper operation (see 'D Clean & store after treatment').
- The adapter, light exit window and the filter of the attachments can become very hot (>210 °C) during usage. Always let the adapter, the light exit window as well as the glass filter and inner parts of the attachments cool down before you touch them.
- This device has a detachable power supply unit. Only use the adapter provided with this device. Do not use an extension cord. The reference number (AD2069x20020HF) can be found on your device. The 'x' in this number refers to the plug type for your country.
- Never leave the device unattended when it is plugged in.
- Always unplug the product when it is fully charged.
- Do not expose the device to direct sunlight or to high temperatures (near hot stoves, in microwave ovens or on induction cookers). Batteries may explode if overheated.
- If the battery is damaged or leaking, avoid contact with skin or eyes. If this occurs, immediately rinse well with water and seek medical attention.
- If the device behaves in a different way during your normal routine (e.g. abnormal smell or vibrations, longer charging time) stop using and charging your device and contact customer support for help.
- Follow the recommended treatment schedule. Increasing the frequency, especially in the initial phase, can increase the risk of skin reactions.


How IPL works

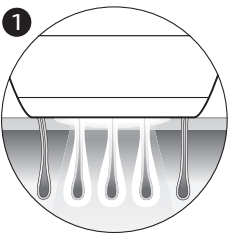With Intense Pulsed Light technology, gentle pulses of light are applied to the skin and absorbed by the hair root. The lighter the skin and the darker the hair, the better the pulses of light are absorbed.

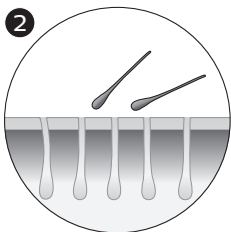The pulses of light stimulate the hair follicle to go into a resting phase. As a consequence, the hair sheds naturally and hair growth is prevented.

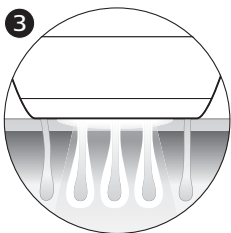The cycle of hair growth consists of different phases. IPL technology is only effective when the hair is in its growing phase. Not all hairs are in the growing phase at the same time. This is why we recommend you to follow the initial treatment phase (4 treatments, every treatment 2 weeks apart) and then the follow-up treatment phase (touch-ups every 4 weeks) to make sure all hairs are effectively treated in the growing phase.

Tip: To assure long lasting hair removal, touch-ups every 4 weeks are recommended.

Suitable body hair colors

Treatment with Lumea is not effective if you have light blond, grey, red or white hairs as hairs with these colors do not absorb enough light.

Note: To check if your body hair color allows usage of the device, consult the hair color table on the foldout page.

Recommended treatment schedule

Initial phase

To achieve effective hair reduction, follow the initial treatment schedule.

Initial treatment schedule


For the first 4 treatments, use Lumea once every 2 weeks to ensure that all hairs are treated.

- Perform a treatment within 3 days before or after the planned treatment date.
- If you have not treated yourself within 3 days before or after the planned treatment date, restart the initial phase. Plan in 4 treatments, 1 treatment every 2 weeks.

Note: If you do not restart the initial phase, you may not see effective hair reduction.

Note: Replacing one of the IPL treatments with another hair removal method (waxing, epilating, etc.) will not help to reach the desired hair reduction.

If you want to remove hairs in between the Lumea treatments, you can use your normal hair removal method (excluding hair removal creams).

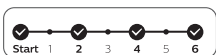Touch-up phase

After the initial phase (4 treatments), we recommend touch-ups every 4 weeks.

Touch-up treatment schedule

Perform a treatment every 4 weeks. Repeat this 8 times to achieve effective hair reduction. This is to maintain results and enjoy smooth skin for months.

- We advise you to perform a treatment within 4 days before or after the planned treatment date.
- The results may vary based on your individual hair growth and also across different body areas.
- If you have completed the touch-up phase and you see the need to continue treatment (hair starts coming back), we advise you to restart the initial phase.

Note: Using the device more often does not increase efficacy.

Tip: Mark the treatment schedule on your calendar to remind yourself of the planned treatments.

What to expect

After initial phase


- After the first treatment, it can take 1 to 2 weeks for the hairs to fall out. In the first weeks following the initial treatments, you still see some hairs growing. These are likely to be hairs that were not in their growing phase during the first treatments.
- After 2-3 treatments, you should see a noticeable reduction in hair growth. However, to effectively treat all hairs, it is important to keep on treating according to the recommended treatment schedule.
- After 4 treatments, you should see a significant reduction of hair growth in the areas that you treated with Lumea. A reduction of hair density should be visible as well.

During touch-up phase


- Keep on treating with frequent touch-ups (every 4 weeks) to maintain the result.

How to use your Lumea Series 9900

A Remove hair & clean your skin

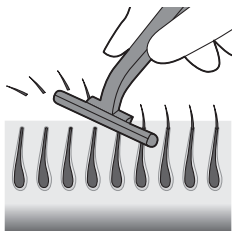
- 1 Remove hair before every treatment as long as hair is still visible. You can either shave, epilate or wax. Make sure you start to treat before new hair is visible on your skin. If you choose to wax, wait 24 hours before you use Lumea.

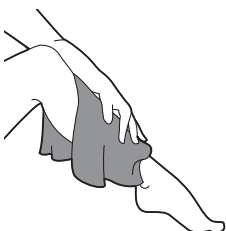- 2 Clean and dry your skin. Make sure that it is free from any lotion or gel.


B Test skin 24 hours before treatment

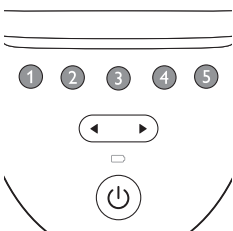
- 1 Test your skin to find the correct light setting:

- For every new body area.
- After recent tanning.

- 2 Select the correct attachment to test your skin (see 'Attachments').

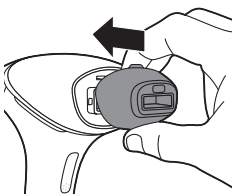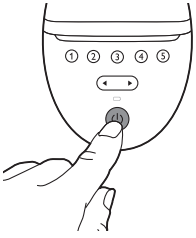
3 Press the on/off button  to turn on the device.

4 Choose an area for your skin test.

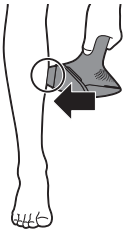
5 Press the device firmly onto your skin and hold it in place.

The integrated safety system prevents unintentional flashing without full skin contact.

6 Confirm that the 'READY to flash' light lights up white to ensure that there is full contact with your skin.

- If the 'READY to flash' light is orange, your skin is too dark to treat. Try Lumea on a lighter body area.

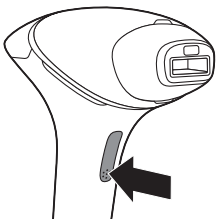
7 Try each setting on your skin and find your comfortable setting.

If the flash feels uncomfortable and you experience pain, lower the setting.

- When the 'READY to flash' light is white, press the flash button on the handle once, to treat one spot with setting ①.


The device makes a soft popping sound. The flash gives you a warm sensation.

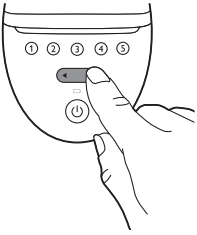- You will see a flash of light.

Note: You can toggle through the light intensity settings by pressing the toggle button.

- Move to next area. Treat one spot with setting ②
- Move to next area. Treat one spot with setting ③
- Move to next area. Treat one spot with setting ④
- Move to next area. Treat one spot with setting ⑤

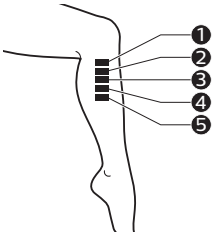24 hours

8 Wait 24 hours.

Note: Even though you found your comfortable setting during the skin test, it may occur that skin reactions show within 24 hours after usage.

9 Check your skin for any reactions (e.g. redness, irritation, rash).


If there are no skin reactions after 24 hours, start the treatment with the highest comfortable setting.

If the highest comfortable setting resulted in a skin reaction, choose a lower setting that did not result in any skin reactions.


If you did not test a lower setting, repeat the skin test.

C Treatment


IMPORTANT: First finish sections A (Remove hair & clean your skin) and B (Test skin 24 hours before treatment) before starting treatment.

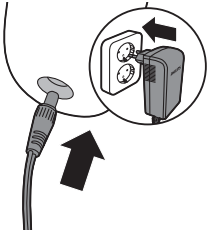
- 1 Connect the power adapter to the device and plug into a power outlet or charge the device and use cordless.

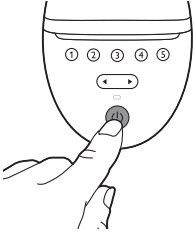
- 2 Press the on/off button  to turn on the device.

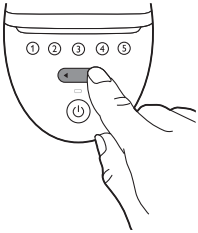- 3 Choose the setting you selected after the skin test (section B, step 9).

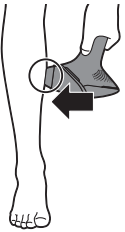
- 4 Press the device firmly onto your skin and hold it in place. The 'READY to flash' light should light up white.

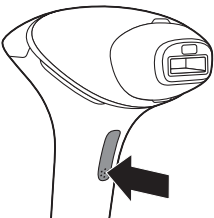- 5 Press the flash button to release a flash on your skin. The flash treats the skin in that spot.

The device makes a soft popping sound. The flash gives you a warm sensation.

After each flash, it takes up to 3.5 seconds until the device is ready to flash again. You can release a flash when the 'READY to flash' light lights up white.

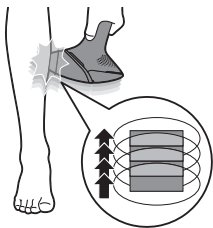
- 6 Move to an adjacent spot and press the flash button again. The surface of the attachment that is in contact with your skin is slightly larger than the light exit window. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area.

DO NOT flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase efficacy.


- 7 Continue until you have treated the full area. There are two treatment modes (see 'Two treatment modes: Stamp & Flash and Slide & Flash').
- 8 Turn off the device when you finished the treatment.
- 9 Remove hair before every next treatment, as long as hair is still visible.

Note: Do not try the device on difficult or sensitive areas (ankle and bony area).


D Clean & store after treatment

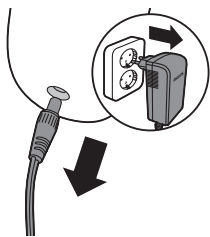
- 1 Unplug the device from the power outlet if you have used it with the power adapter.
- 2 Let the device cool down.

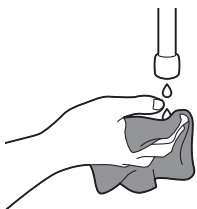
- 3 Moisten the soft cloth supplied with a few drops of water.

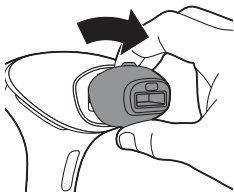- 4 Take the attachment off the device.

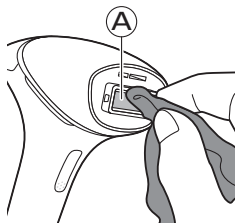- 5 Clean the light exit window of the device (A).

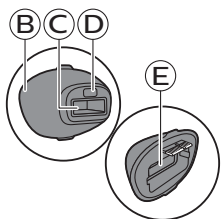- 6 Clean all surfaces of the attachment, including the metal reflector and window of the attachment (parts B, C, D, E).
- 7 Let all parts dry thoroughly.
- 8 Store in a dust-free place at room temperature.

Learning more about your Lumea Series 9900

Charging

Fully charge the batteries before you use the device for the first time and when they are empty. Charging the batteries takes up to 1 hour and 40 minutes. Charge the device when the charging light lights up orange during use to indicate that the battery is low and will run out soon.

Fully charged batteries provide at least 130 flashes at light intensity 5.

Charging the device

- 1 Turn off the device.
- 2 Insert the small plug into the device and the adapter into the power outlet.
 - During charging, the charging light flashes white.
 - When the batteries are fully charged, the charging light lights up white continuously.
 - Never cover the adapter and device during charging.

Note: The adapter and device can feel warm during charging. This is normal.

Note: This device is equipped with battery-overheat protection and does not charge if the room temperature exceeds 40 °C.

- 3 After charging, remove the adapter from the power outlet and pull the small plug out of the device.

Tip: Charge the device after each use to save battery life.


Note: The battery capacity is not enough for a full body treatment. We advise you to use the device corded while treating large body areas such as legs.

Attachments


For optimal results, it is important to use the correct attachment intended for the specific body area you want to treat. Lumea offers full-body treatment with several attachments designed for use on specific body areas. **SenseIQ** technology adapts the treatment program when a different attachment is connected.

Note: The device may not work anymore and show an error when there is dirt on the attachment connector. Clean the contact leads when this occurs.

To place the attachment, simply snap it onto the light exit window.

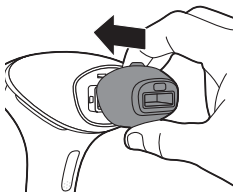To remove the attachment, pull it off the light exit window.

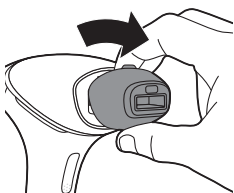Body attachment

The body attachment has the largest treatment window and a curved-in design to effectively cover and treat areas below the neckline, especially large areas such as legs, arms and stomach.

Facial attachment

The facial attachment has a precise flat design with extra integrated red glass filter for safe and precise treatment on the sensitive skin above the upper lip, chin and sideburns. The device is only intended for removing unwanted body hair from areas below the cheekbones. When you use the device on areas with sharp contours like the jawline and chin, it may be difficult to establish full skin contact and release a flash.


Tip: Place your tongue between your upper lip and your teeth or puff out your cheeks to make the treatment easier.

Underarm attachment (BRI976, BRI977)

The underarm attachment has a specifically curved-out design to treat underarm hairs that might be hard to reach.

Bikini attachment (BRI976, BRI977)

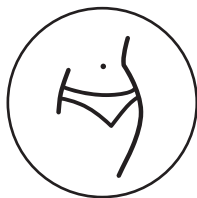The bikini attachment has a specialized design for effective treatment of the bikini area. It has a curved-out design and has a higher light intensity compared to body attachment. Hair in this area tends to be thicker and stronger than leg hair.

Precision attachment (BRI973)

The precision attachment has a curved-out design for use on bikini and underarm areas. It has a medium sized window with a glass window. It is designed for precise and effective coverage of bikini and underarm areas.

SmartSkin sensor

Lumea offers **SenseIQ** technology for a personalized hair removal treatment. The integrated SmartSkin sensor measures the skin tone at the start of each session and during the session. This provides two features:

- 1 The SmartSkin sensor prevents you to treat areas of your skin that are too dark: if it detects a skin tone that is too dark for treatment with Lumea, the 'READY to flash' light starts blinking orange. It does not flash when you press the flash button. The device automatically disables. This prevents you from developing skin reactions.
- 2 The SmartSkin sensor will help you to select the light intensity setting that is comfortable for your skin. For this you can use the SmartSkin button .

Selecting the light intensity

Lumea provides 5 different light intensities, you will always be able to change to a light intensity setting you find comfortable.

- 1 Press the on/off button  to turn on the device. The device starts to operate at light intensity setting 1.

Note: To adjust the light intensity setting manually without using the SmartSkin button, use the toggle button until you have reached the required setting. The corresponding intensity light level lights up white.

- 2 Press the SmartSkin button  and put the device to your skin.
- 3 After the device has scanned your skin, the light intensity indicators will blink white to indicate which intensities are comfortable based on your skin tone.
- 4 Press the confirmation button  to confirm that you want to use the indicated setting.

Lumea automatically uses the highest indicated intensity, which is shown by a blinking intensity indicator.

Lumea gives you the freedom to select the light intensity that you find comfortable. Skin sensitivity can differ per person. Therefore, the skin test is the most important guide for selecting your light intensity setting (see 'B Test skin 24 hours before treatment').

- 5 Using Lumea should never be uncomfortably painful. If you experience discomfort, reduce the light intensity setting. You can do this by using the toggle buttons.


Note: The device is automatically disabled when your skin tone is too dark (dark, according to the skin tone table), to prevent you from developing skin reactions. The 'READY to flash' light will blink orange to indicate if the skin tone is too dark.

- 6 After moving to another body area or after recent tanning (see 'Using your Lumea Series 9900 before and after tanning'), perform a skin test to determine a comfortable light intensity setting. To re-activate the setting indication feature, press the SmartSkin button .

Note: Perform a skin test before you start to treat another body area or after a skin tone change. Determine the light intensity setting for each body area separately.


Two treatment modes: Stamp & Flash and Slide & Flash

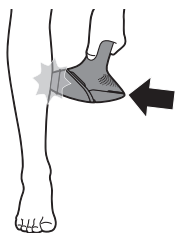Your Philips Lumea has two treatment modes for more convenient use on different body areas:

- The Stamp & Flash mode is ideal to treat small or curvy areas like knees and underarms. Simply press and release the flash button to release a single flash.

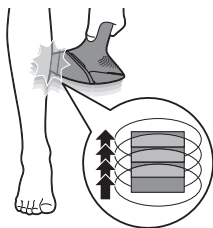- The Slide & Flash mode offers convenient use on larger areas like legs. Keep the flash button pressed while you slide the device over your skin to release several flashes in a row.

Using your Lumea Series 9900 before and after tanning

Tanning with natural or artificial sunlight

Intentionally exposing your skin to natural or artificial sunlight with the aim of developing a tan influences the sensitivity and color of your skin. Therefore the following is important:

- After each treatment, wait at least 48 hours before tanning. Even after 48 hours, make sure that the treated skin does not show any redness from the treatment anymore.
- If you expose your skin to the sun (without tanning intentionally) in the 48 hours after treatment, use a sunblock SPF 50 on the treated areas. After this period, you can use a sunblock SPF 50 for two weeks.
- After tanning, wait at least 2 weeks before you use Lumea and perform a skin test. This is to determine the appropriate light intensity setting (see 'B Test skin 24 hours before treatment').
- Do not use Lumea on sunburned body areas.

Note: Occasional and indirect sun exposure does not qualify as tanning.

Tanning with creams

If you have used an artificial tanning lotion, wait until the artificial tan is completely gone before you use the device.

Traveling with the device

When you travel with the device, check with the airplane carrier to confirm that the device can be carried and/or used on the airplane.

After use

Aftercare

After use, you can apply lotions, creams, deodorant, moisturizer or cosmetics to the treated areas.

Warning: If you experience skin irritation or skin redness after treatment, wait until it disappears before applying any product to your skin. If you experience skin irritation after applying a product to your skin, wash it off with water. If you still get skin reactions, stop using the device and consult your doctor.

Common skin reactions

Your skin may show slight redness and/or may prickle, tingle or feel warm. This reaction disappears quickly.

Dry skin and itching may occur due to shaving or a combination of shaving and light treatment. You can cool the area with an ice pack or a wet face cloth. If dryness persists, you can apply a non-scented moisturizer on the treated area.

Rare side effects

- Burns, excessive redness (e.g. around hair follicles) and swelling: these reactions occur rarely. They are the result of using a light intensity that is too high for your skin tone. If these reactions do not disappear within 3 days, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Skin discoloration: this occurs very rarely. Skin discoloration manifests itself as either a darker patch (hyperpigmentation) or a lighter patch (hypopigmentation) than the surrounding area. This is the result of using a light intensity that is too high for your skin tone. If the discoloration does not disappear within 2 weeks, we advise you to consult a doctor. Do not treat discolored areas until the discoloration has disappeared and your skin has regained its normal skin tone.
- Skin infection is very rare but is a possible risk following a (micro)wound, a skin burn, skin irritation, etc.
- Epidermal heating (a sharply defined brownish area which often occurs with darker skin tones and is not accompanied with skin dryness): this reaction occurs very rarely. In case this reaction does not disappear within 1 week, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.
- Blistering (looks like small bubbles on the surface of the skin): this occurs very rarely. In case this reaction does not disappear within 1 month or when the skin gets infected, we advise you to consult a doctor. Wait with the next treatment until the skin has healed completely and make sure you use a lower light intensity.

- Scarring: often the secondary effect of a burn, which can take longer than a month to heal.
- Folliculitis (swelling around hair follicles combined with pustule formation): this reaction occurs very rarely and is the result of bacteria penetrating the damaged skin. In case this reaction occurs, we advise you to consult a doctor as folliculitis may need antibiotic ointment.
- Pain: this can occur during or after treatment if you use the device on skin that is not hair-free, if you use the device at a light intensity that is too high for your skin tone, if you flash the same area more than once and if you use the device on open wounds, inflammations, infections, tattoos, burns, etc.

Troubleshooting

Problem	Possible cause	Solution
The device/adapter becomes warm during use.	It is normal for the device and adapter to become warm (but not too hot to touch) during use.	Use the device in a slightly cooler environment and/ or let it cool down before continuing use.
	The fan is not working.	Check if the attachment is well connected. Clean the connectors on the attachment if needed. In case the attachment is well connected and it is possible to produce flashes with the device, contact the Consumer Care Center in your country, your Philips dealer or a Philips service center.
	The cooling airflow of the fan is blocked by hands or a towel.	Make sure the cooling airflow of the fan is unblocked.
When I place the device on the skin, it does not release a flash. The 'READY to flash' light blinks orange.	Your skin tone in the area to be treated is too dark.	Treat other body areas with lighter skin tones with Lumea.
The 'READY to flash' light blinks orange and all 5 intensity lights blink as well.	The device needs to be reset.	To reset the device, take the plug out of the socket, wait for 30 minutes to let the device cool down. The device should function normally again. In case it does not work again, contact the Consumer Care Center in your country.
The 'READY to flash' light does not light up white.	The device is not completely in contact with your skin.	Place the device at a 90 degree angle on the skin so that the integrated safety system is in contact with your skin.

Problem	Possible cause	Solution
The device produces a strange smell.	The light exit window or the SmartSkin sensor is dirty.	Clean the light exit window and the SmartSkin sensor carefully.
	You have not removed the hairs on the area to be treated properly. These hairs may get burned and can cause the smell.	Pretreat your skin before you use Lumea.
The device does not flash, the fan is not switched on and all 5 intensity lights blink.	The attachment is not attached properly.	Make sure you attach the attachment completely. If necessary, clean the electronic contacts on the attachment.
The skin feels more sensitive than usual during treatment. I experience discomfort when I use the device.	The light intensity setting you use is too high.	Check if you have selected a comfortable light intensity setting. If necessary, select a lower setting.
	You did not remove the hairs on the areas to be treated.	Pretreat your skin before you use Lumea.
	The UV filter of the light exit window is broken.	If the UV filter is broken, do not use the device anymore. Contact the Consumer Care Center in your country, your Philips dealer or a Philips service center.
You treated an area for which the device is not intended.	<p>Never use the device on irritated (red or cut), sunburned, recently tanned or fake-tanned skin. Never use the device if you have infections, eczema, burns, inflammation of hair follicles, open lacerations, abrasions, herpes simplex (cold sores), wounds or lesions and haematomas (bruises) in the areas to be treated.</p> <p>Never use the device on the following areas: on moles, freckles, large veins, darker pigmented areas, scars, skin anomalies without consulting your doctor. Never use the device on warts, tattoos or permanent make-up, around the eyes and on or near the eyebrows, on the lips, perineum, scalp, nipples, areolas, labia minora, vagina, anus and the inside of the nostrils and ears. Men must not use it on the face and neck including all beard growing areas and whole genital area. This can result in a burn and a change in skin color, which makes it potentially harder to identify skin-related diseases.</p>	
There is no glass window in my attachment.	This is normal.	No action required: there is no glass window in the body and underarm attachments. Only the facial attachment has a red glass filter. The precision and bikini attachments have a transparent glass window.

Problem	Possible cause	Solution
The skin reaction after treatment lasts longer than usual.	You have used a light intensity setting which is too high for you.	Select a lower intensity next time (see 'Selecting the light intensity').
The flash is very bright to my eyes. Do I need to wear goggles?	No, Philips Lumea does not hurt your eyes.	The scattered light produced by the device is harmless to your eyes. Do not look at the flash while using the device. It is not necessary to wear goggles during use. Use the device in a well-lit room so that the light is less glaring to your eyes. Be sure to make good skin contact to avoid scattered light.
The results of the treatment are not satisfactory.	You have used a light intensity setting which is too low for you.	Select a higher setting next time.
	You did not cover all areas to be treated.	The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase efficacy.
	The device is not effective on your body hair color.	If you have light blond, grey, red or white hair, the treatment is not effective.
	You do not use the device as often as recommended.	To remove all hairs successfully, we advise you to follow the recommended treatment schedule. You can reduce the time between treatments, but do not treat more often than once every two weeks.
	You respond more slowly to IPL treatment.	Continue using the device for at least 6 months, as hair growth can still decrease over the course of this period.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU).

- This symbol means that this product contains a built-in rechargeable battery which shall not be disposed of with normal household waste (2006/66/EC). Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery. Philips Lumea products contain recyclable materials and should not be put into the municipal waste stream. Refer to the Philips Lumea support website for recycling options. Do not dispose of in fire.

Technical specifications

Transport and storage conditions	Lumea remains operational within its specification in NORMAL USE after transport or storage in the following environmental conditions: -25 °C to 70 °C up to 90 % RH, non-condensing, 700 hPa to 1060 hPa
Operating conditions	Lumea complies with its specifications when operated in NORMAL USE under the following environmental operating conditions: 5 °C to 35 °C up to 90 % RH, non-condensing, 700 hPa to 1060 hPa
Emitted wavelengths	565 to 1400 nm
Pulse duration	1–2.5 milliseconds, for all settings
Flash interval	Between 1–3.5 seconds depending on battery/mains connected use
Optical exposure	Body attachment: 2.5–5.6 J/cm ² Underarm attachment: 2.5–6.5 J/cm ² Precision/Bikini attachment: 2.5–6.5 J/cm ² Facial attachment: 2.4–5.9 J/cm ²
Adapter input	100–240 V ~ 50/60 Hz 1.5 A
Adapter output	19.5 V === 3.33 A 65.0 W
Average active efficiency	≥88.0 %
Efficiency at low load (10 %)	≥79.0 %
No-load power consumption	<0.21 W
Adapter type number	The reference number (AD2069x20020HF) can be found on your device. The letter 'x' denotes type of plug, it can be: 0, 1, 2, 3, 5, 7, 9, or H.

Electromagnetic compatibility - Compliance information

This device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Explanation of symbols

The warning signs and symbols are essential to ensure that you use this product safely and correctly and to protect you and others from injury. Below you find the meaning of the warning signs and symbols on the label and in the user manual.

This symbol identifies the manufacturer of the product: Philips Consumer Lifestyle B.V., Tussendiepen 4, 9206 AD Drachten, The Netherlands.

This symbol means that you have to read the user manual carefully before you use the device. Save the user manual for future reference.

This symbol means that the user should consult the instructions for use for important safety information such as cautions that cannot, for a variety of reasons, be presented on the device itself.

This symbol means that the user should consult the instructions for use for important safety information such as contraindications and warnings that cannot, for a variety of reasons, be presented on the device itself.

This symbol indicates the manufacturer's serial number.

This symbol identifies the manufacturer's catalogue number.

This symbol identifies that the equipment meets the safety requirements specified for Class II equipment.

This symbol indicates that the equipment is suitable for direct current only.

This symbol indicates that the device should not be used around the eyes and on or near the eyebrows.

This symbol means: 'Do not use this device in a bathtub, shower or water-filled reservoir.'

This symbol means that the transformer is designed to be short-circuit proof.

This symbol means that the transformer is considered a Switch Mode Power Supply (SMPS) unit.

This symbol means that the device is designed for indoor use.

t_a

This symbol indicates the rated maximum ambient temperature at which the transformer may be operated continuously under normal conditions of use.

≤ 3000m

This symbol indicates that the appliance is intended to be usable up to the maximum altitude 3000 m.

A recycling symbol to identify the material of which an item is made, to facilitate recycling or other reprocessing. The symbol may include a number and/or abbreviation.

This symbol indicates that a specific separate supply unit is required for connecting electrical equipment to the supply mains.

This symbol identifies the switch or switch position by means of which part of the equipment is switched on in order to bring it into the stand-by condition, and to identify the control to shift to or to indicate the state of low power consumption.

IP22

This symbol indicates the protection against ingress of solid foreign objects ≥ 12.5 mm and against harmful effects due to the ingress of vertical dripping water when tilted at 15 degrees.

IP30

This symbol indicates the protection against ingress of solid foreign objects ≥ 2.5 mm and that there is no protection for harmful effects due to the ingress of water.

This symbol indicates the energy efficiency level.

This symbol indicates that the adapter should not be used if the plug pins are damaged.

This symbol means an object is capable of being recycled - not that the object has been recycled or will be accepted in all recycling collection systems.

This symbol indicates WEEE, waste electrical and electronic equipment. Electrical waste products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice (see 'Recycling').

This symbol means that this product contains batteries which shall not be disposed of with normal household waste (2006/66/EC).

The symbol for the Green Dot ('Der Grüne Punkt' in German). The Green Dot is the license symbol of a European network of industry-funded systems for recycling the packaging materials of consumer goods.

Symbol applied to products to indicate that they conform with relevant EU directives regarding health and safety or environmental protection.

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