CURAPROX HYDROSONIC BLACK IS WHITE

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Congratulations

on purchasing a Hydrosonic «Black Is White»
sonic toothbrush.

«Black Is White» from CURAPROX: oral health with style.

Whitening bristles: thanks to the activated carbon
particles in the bristles, the «Carbon» brush head
effectively removes discolouration without damaging
the enamel or gums.

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English
see page 18

Read this user manual
carefully before you use the
appliance and save it for
future reference

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Warning

If you've recently had surgery on your teeth or gums, consult your dentist before using the Hydrosonic «Black is White». If you experience excessive bleeding after using this device, please contact your dentist immediately. This toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance. The Hydrosonic «Black is White» complies with safety standards for electromagnetic devices. If you wear an implant or pacemaker and have questions about the device, please contact your doctor or the manufacturer of your implant or pacemaker.

Safety note

Warning: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. This appliance contains batteries that are non-replaceable. Before using, please check that the voltage

stated on the device complies with your local mains voltage. If the Hydrosonic charger or connecting cable is damaged, please do not use or attempt to repair it yourself. The power cord cannot be replaced. If the cord is damaged, the entire device must be replaced. All repairs must be handled by a dealer or our customer service. Keep the charger away from water, and do not immerse in water or other liquids. Cleaning the charger: Disconnect from the mains before cleaning and reconnect only once the charger is completely dry. Never charge the device unattended.

Environmentally friendly disposal

Do not dispose of this product with normal household waste at the end of its service life. Take it to a recycling centre or official collection point. This helps protect the environment. The integrated rechargeable battery contains substances that may endanger the environment. So please take the used handset to an official recycling centre or collection point.

Proper disposal of this product

Within the EU, this symbol indicates that this product must not be disposed of with household waste. Old appliances contain valuable recyclable materials. They should be recycled in order to prevent harm to the environment or human health through improper waste disposal. Therefore, please dispose of old devices through a proper collection system, or send the device back to where you purchased it for proper disposal. They will then ensure that the device is recycled accordingly.

Technical data for USB plug

Adapter Input: 100-240 VAC, 50/60Hz, 0.2A max.

Adapter Output: 5 VDC, 600mA

Technical data for charging station with USB cable

Input: 5 VDC, 600mA

Single base, rubber feet

Technical data for handset

- Battery: lithium-ion (3.7 VDC, 700mAh)
- Optimal battery charging time: 10 hours
- Operating time per battery charge: approx. 60 mins.
This is equivalent to 15 days if you brush your teeth for four minutes each day.
- Brushing movements: 22,000, 32,000 to 42,000
- Quadrant timer: 30 seconds
- Brushing timer: two minutes
- Automatic switch-off: after four minutes
- Accessories: 2× «Carbon» brush head
 - 1× charger with USB cable
 - 1× USB plug
 - 1× travel case
 - 1× instruction manual

Conformity

Troubleshooting

Problem	Possible cause	Solution
• The device isn't working	• The device is not charged	• Place the handset securely on the charger • Charge the Hydrosonic completely before first use. The charging takes about 10 hours.
• Brushing performance is getting weaker	• The brush head is dirty or worn	• Clean or change the brush head
• You feel a slight tingling when using the sonic toothbrush	• If you are using a sonic toothbrush for the first time, this is a normal sensation	• You will become accustomed to the sensation after using the sonic toothbrush for a while

HYDROSONIC BLACK IS WHITE


Getting started	20
Charging: extra long life	20
«start», «clean» and «smile» modes	21
Movement: up and down	22
Jaw quadrants: always the same	23
Whitening bristles	24
Ergonomics: CURACURVE®	25
Oral health with style	26
Customer service	27

Register now
and get an extra
year guarantee.

[www.curaprox.com/
warranty](http://www.curaprox.com/warranty)

Contents


- 1× Hydrosonic «Black Is White»
- 2× «Carbon» brush head
- 1× charger with USB cable
- 1× USB plug
- 1× travel case
- 1× instruction manual

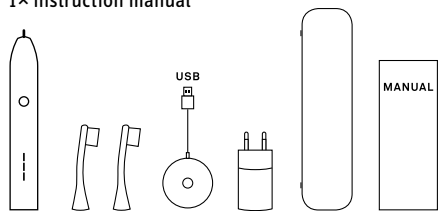GETTING STARTED


Select your brush head

Attach the brush head, moisten with water, add toothpaste. To start and switch off: press the button briefly. Press and hold the button to switch between modes: start, clean or smile.

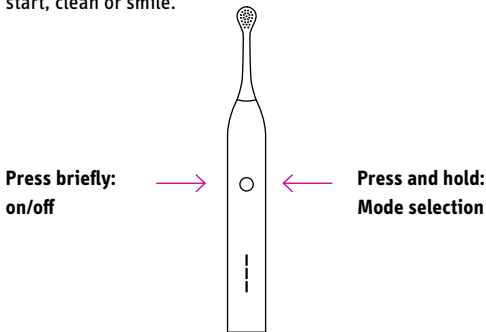Clean without applying pressure: the bristles do the work; simply hold the brush against each tooth for two to three seconds without applying pressure.

After brushing your teeth: rinse and remove the brush head. Or simply leave it on. Caution: do not bend.

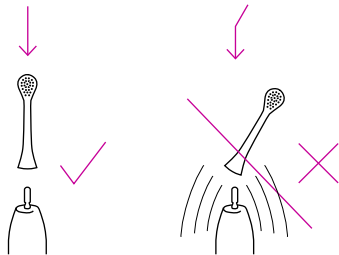CHARGING: EXTRA LONG LIFE

- Charge the Hydrosonic «Black Is White» completely before first use.
- Optimal battery charging time: 10 hours.
- Charging: place the Hydrosonic «Black Is White» into the charger and connect it to the mains. To do this, plug the charger's USB connector into the USB adapter and plug this into the mains. You can also charge the device using the USB port on your laptop or computer.
- Check the charge status: remove the Hydrosonic «Black Is White» from the charger or switch it off.
- Fully charged: all LEDs are lit.

- Half charged: the two bottom LEDs are lit.
- Charged below 20%: the lowest LED lights up red instead of pink.
- When fully charged, the battery lasts a good 60 minutes, which equates to approximately two weeks of brushing for four minutes each day.
- **TIP:** For extra long life: do not charge the device until the battery status falls below 20% (the bottom LED lights up red instead of pink).
- **NOTE:** The cleaning performance is not affected by the charge level. The Hydrosonic «Black Is White» cleans just as well when the battery is low as when it is fully charged.
- **CAUTION:** To clean the device, wipe the handset and charger clean with a damp cloth. Never clean the device with running water.

THE MODES: START, CLEAN, SMILE

When you switch on the Hydrosonic, it will start on the last selected mode.

- The bottom LED is lit up: «start» mode stimulates circulation in the gums, with 22,000 brush strokes per minute. This mode will also help you get used to the slight tingling caused by a sonic toothbrush.
- The middle LED is lit up: «clean» mode with 32,000 brush strokes per minute optimally removes plaque and bacteria.
- The top LED is lit up: «smile» mode is recommended to clean particularly sensitive teeth and gums, with 42,000 brush strokes per minute.

Every mode provides excellent dental cleaning – the «clean» and «smile» modes provide hydrodynamic as well as mechanical cleaning.

MOVEMENT: BACK AND FORTH


Our oral health professionals recommend a simple vertical back-and-forth movement, without applying pressure. This has three advantages: It is gentle on your gumline. It provides particularly good care. It delivers hydrodynamic as well as mechanical cleaning.

Hydrodynamics: turbulent


The bristles move very quickly, generating turbulence in the mixture of saliva, water and toothpaste. This turbulence also reaches areas which cannot be reached by the bristles, thus ensuring that cleaning takes place between your teeth as well.

Effective: slanted

- Position the brush at a 45-degree angle so it lies half on the gums, half on the teeth. This is the ideal way to clean your teeth with the ergonomic hexagonal handle.


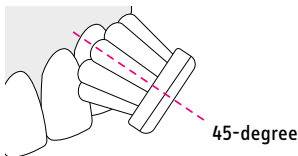
- Simply hold the brush loosely without applying pressure.
- Clean from tooth to tooth: two to three seconds on each tooth.
- While cleaning, slightly tilt the Hydrosonic «Black Is White» once or twice to clean the recesses and spaces between your teeth particularly well.

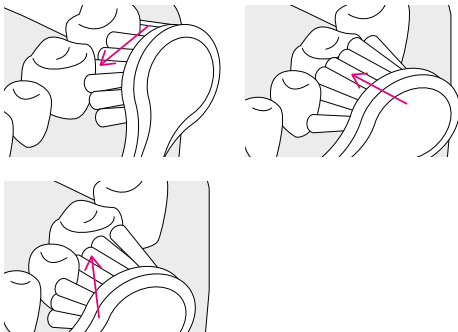JAW QUADRANTS: ALWAYS THE SAME

- After 30 seconds, three short, consecutive pauses of 0.2 seconds each will indicate that it's time to move to a different quadrant in your mouth.

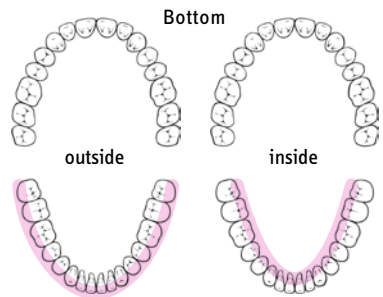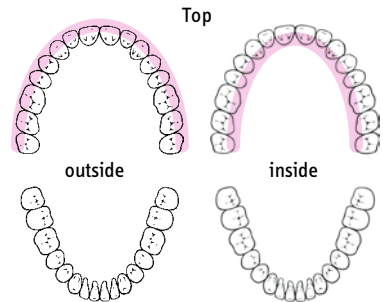- Important: the cleaning process should be repeated for every jaw quadrant – ideally always in the same order. This way, you won't miss out any teeth.
- A long pause indicates that the two minutes have elapsed.
- After four minutes, the Hydrosonic «Black Is White» will switch off.
- You can switch the toothbrush off at any time by pressing the button.

Chewing surfaces: don't forget

Children and adolescents in particular must also clean the chewing surfaces of their teeth due to the often deep fissures.

Dentist recommendation: ideally longer

Dentists recommend a minimum cleaning duration of two minutes. Cleaning a little longer promotes oral health.

WHITENING BRISTLES

The bristles of the «Black Is White» «Carbon» brush head feature soft, highly effective activated carbon particles. These particles allow you to gently remove discolouration, without damaging the enamel.

We have been using CUREN® bristles instead of conventional nylon brushes from the very beginning. The advantage: CUREN® hardly absorbs any water, so performance remains constant. Another particularly important feature: CUREN® is unbelievably soft and feels wonderfully comfortable – even with 42,000 brush strokes per minute.

We are proud of our innovative Swiss brush head technology – and the fact that we produce our brush heads in Switzerland. We are thrilled to be contributing to improved oral health in a simple, comfortable manner.

MADE IN SWITZERLAND

Order replacement brush heads: specialist shops, pharmacies, chemists, www.curaprox.com

ERGONOMICS: CURACURVE®

Yes, our brush heads bend easily. This 10 degree bend is crucial: the ergonomics are significantly improved, difficult areas such as the molars or the inside of the front teeth are easily accessible – the perfect bend thanks to CURACURVE®.

10 degree


0

Now perfectly accessible thanks to CURACURVE®: the inside of your front teeth. Always stay on each tooth for two to three seconds.

CURACURVE® also facilitates cleaning the rear molars at the top and bottom. To do so, guide the Hydrosonic «Black Is White» steeply behind the molar. Yes, for two to three seconds.

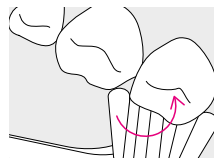We recommend switching the brush head every three months.

ORAL HEALTH WITH STYLE

With our «Black Is White» product line, oral health becomes a lifestyle. The following products are part of the line:

The whitening toothpastes «Black Is White»

- **White teeth:** The toothpastes remove discolouration with the aid of activated carbon without grinding or bleaching.
- **Strong oral health:** Enzymes protect against decay and support the function of saliva. Hydroxyapatite remineralises the enamel and closes exposed dentinal tubules.
- **Plus:** With physical blue filter. With fluoride. No SLS. No triclosan. No bleaching agents. No plastic particles. Usable without restrictions like normal toothpaste. Made in Switzerland..

The Hydrosonic «Black Is White» and the «Black Is White» whitening toothpaste: an unbeatable power-duo.

And for between brushes: «Black Is White» chewing gum

It tastes lemony minty fresh. The coolest part: it contains activated carbon, thus removing discolouration and whitening teeth without damaging them. Hydroxyapatite fills even the smallest beginnings of decay and remineralises enamel. Glucose oxidase enables balanced oral flora and has an antibacterial effect. Xylitol protects against decay and provides a lovely, sweet flavour.

CUSTOMER SERVICE

Guarantee terms

CURADEN AG provides guarantee for material, function and manufacturing defects for two years from the date of purchase. This may be carried out by way of repair work or a replacement device. In order for a guarantee claim to be made, the original purchase receipt must be submitted. Keep your receipt as proof of purchase. The guarantee shall expire if the device is subject to interference or changes to its original state.

Guarantee exclusions

In the following cases, CURADEN AG shall not provide guarantee services: Wear and tear of the device or parts of the device due to normal use, or if the device was not used, not properly maintained, or if the safety instructions provided with the device were not observed. Repairs may only be carried out by service partners authorised by CURADEN AG.

Service

In order to make a guarantee claim, send the device and original receipt to your country's official CURAPROX representative. If you do not know who your country's representative is, visit our website www.curaprox.com or contact us directly in Switzerland: ask@curaprox.com

- Please discuss any questions about dental or orthodontic care with your dentist. You can also contact us, especially if you have technical problems or if something is not working: ask@curaprox.com
- Our FAQs are often helpful: www.curaprox.com/faq Here, you will also find out how to exchange your Hydrosonic «Black is White» if necessary.

CURAPROX «Black Is White» products are available in specialist shops, pharmacies and chemists as well as at www.curaprox.com

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