

BLACK WHITENING®

GEBRAUCHSANWEISUNG

MANUEL D'UTILISATION

INSTRUCTION MANUAL

megasmile

ENGLISH



FOR SAFETY REASONS

Warning:

- If you have had oral or gum surgery lately, please consult your dentist. Consult your dentist if excessive bleeding occurs after using this tooth brush or if bleeding continues.
- Do not take the device apart, (Exception: Disposal).
- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe manner, and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance.
- Keep charger away from water or other liquids! Do not leave it near a bathtub or a basin.
- Keep charger away from hot surfaces and direct sunlight if connected to the socket.


- The Black Whitening® toothbrush complies with safety standards for electromagnetic devices.


If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device.

Please read this instruction manual carefully before using the device and save it for future reference.

THE TECHNOLOGY OF BLACK WHITENING®

Black Whitening® is a rechargeable hydro-active toothbrush with sonic-technology for long-lasting healthy and dazzling white teeth. Surges of the sonic technology create pulsating flows of fluids, which easily enter even the narrowest dental spaces. The microparticles of charcoal efficiently remove stains and polish tooth surfaces without harming the enamel.



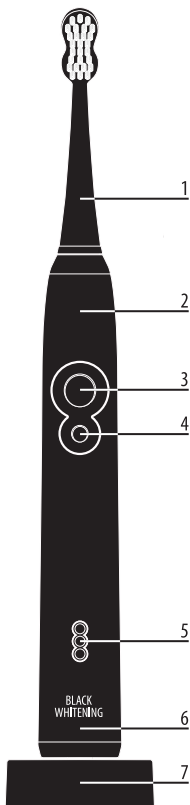

OVERVIEW

- 1| Sonic brush head
- 2| Slip-safe handle with brushing timer and quadrant timer
- 3| On/off button
- 4| Brushing modes with memory (Mode Selection):
SENSITIVE - INTENSE - DEEP CLEANSING mode
TRAVEL LOCK (press for 4 sec.)
- 5| LED display: Indicates the selected clean mode and state of charge.
- 6| Rechargeable Li-Ion battery
- 7| USB Charger with automatic voltage adjustment 100-240V

OPERATION

Align the brush head so the bristles face the front of the handle. Firmly press the brush head down on the metal shaft until it clicks into place.

Use the USB cable and connect the charger with the external power adapter. Plug the adapter into a socket and place the hand piece onto the charger. The LED display will light up for about 10 sec. Afterwards all LED-lights will switch off. The Black Whitening® toothbrush is charging.

Pull the brush head off the handle for replacement. Don't twist the brush head when removing it from the handle as it may damage the motor. Don't immerse the toothbrush handle under water and do not leave the handle on wet surfaces. Dry the toothbrush handle off with a soft cloth after brushing.

CHARGE - LEVEL INDICATOR

The charging takes about 12 hours. Charge the Sonic Generation® toothbrush for 24 hours before first use.

The LED display will blink for 10 sec. if the device is removed from or placed on the charger, or switched off. Three blinking LED-lights indicate a fully charged battery. If one single light is blinking the battery is charged less than 25%.

Don't keep the toothbrush handle on the charger base all the time. Charge the toothbrush only when the charge-level is low, to maximise battery life.

USAGE OF THE BLACK WHITENING®

Wet the bristles and apply a small amount of toothpaste. Gently place the brush against the teeth at a slight angle towards the gum line.

The toothbrush remembers the last chosen mode and will start automatically by pressing the on button in this mode (Memory). Repeatedly pushing on the button MODE SELECTION will change the clean mode:





SENSITIVE Mode:

Upper LED-light glows



INTENSIVE Mode:

Upper two LED-lights glow



DEEP CLEANSING Mode:


All three LED-lights glow

Move the toothbrush in an up-and-down manner from tooth to tooth, without exertion; let the brush do the work. Rest 2-3 sec. on one tooth before moving to the next one. Make sure that the longer bristles penetrate the dental spaces between teeth. For each quadrant, move the brush first along the internal tooth side over to the chewing surface, and at last along the external tooth side.



Rest 2-3 sec. on one tooth before moving to the next one.

After 30 sec., the toothbrush signals to change the quadrant by a short interruption in its vibration. Repeat this cleaning process for each quadrant. When the time has run out the device automatically deactivates. If needed, it is possible to pursue the brushing by simply reactivating the device.

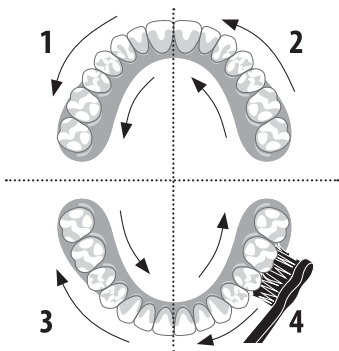

Important:

do not exert pressure on the toothbrush. Gentle surface contact along the tooth alignments will lead to better cleaning results and maintain the health of your teeth and your gum line.

Your Black Whitening® is also safe to use on: braces, implants and dental restorations (fillings, crowns, veneers).

Replace Black Whitening® brush heads every 3 months for hygienic reasons. Brush heads wear out sooner when used on braces.

Ask your dentist, dental hygienist or assistant after every appointment if and how you are supposed to use your sonic toothbrush. Dental professionals are always competent for follow-up questions.

CLEANING MODES

The Black Whitening® is made for different uses. Depending on the user's needs, it is possible to change the clean mode.

Select the different modes by repeatedly pushing the MODE SELECTION button.

The Black Whitening® remembers the last chosen mode so that it does not have to be selected the next time you activate the device (Memory function).

SENSITIVE mode

By selecting SENSITIVE mode the upper LED-light glows. The sonic waves that move with a speed of up to 45.000 strokes per minute allow the bristles to vibrate gently. This is a gentle way of cleaning enamel and exposed dental necks. We recommend the SENSITIVE mode to people with sensitive teeth, exposed dental necks, periodontitis, implants, bridges, and other dental replacement solutions. This mode is ideal for regular care of teeth with fluoride gel.

INTENSIVE mode

By selecting INTENSIVE mode the two upper LED-lights glow. The surges of the sonic technology create a pulsating flow of fluid which easily enters even tiny dental spaces, neutralizes debris, and assures a perfect cleaning all around teeth.

DEEP CLEANSING mode

By selecting DEEP CLEANSING mode all three LED-lights glow. Increased brushing movements on the tooth surfaces allow increased cleaning. Teeth are freed of unpleasant discoloration. We recommend the DEEP CLEANSING -mode especially in case of discoloration as a result of nicotine, wine, tea or coffee, and for occasional extra-thorough teeth cleaning.

Quadrant timer

The quadrant timer is a 30 second interval timer that briefly interrupts the brushing to signal the recommended change to the next quadrant. It reminds you to brush the four-quadrants of your mouth equally and thoroughly, for a better overall-cleaning.

Brushing timer

Dentists recommend a minimum of 2 minutes brushing time. After the recommended brushing time, depending on when the selected mode; has been reached, the Black Whitening® deactivates.

TRAVEL LOCK

The Black Whitening® - travel lock; locks and unlocks the buttons in order to prevent accidental operation during travel.

Locking the device:

In order to lock the buttons, push the button: MODE SELECTION for about 4 sec., until the LED-lights blink. The device is now safe from accidental activating.

Unlock:

To unlock the device push the MODE SELECTION button for about 4 sec. The LED-lights show by blinking that the device is unlocked. An alternative is to put the toothbrush back on the charger base. This unlocks it automatically.

QUESTIONS AND ANSWERS

I have sensitive teeth and gums.

What is the best way to brush with my Black Whitening®?

If you have sensitive teeth and gums, plan for a 1-2 week adjustment period to your Black Whitening®. During this time, we recommend using the SENSITIVE mode.

How can I best remove discolorations?

We recommend the DEEP CLEANSING mode in order to rid teeth of discolorations.

What do I have to consider when I brush my teeth with Black Whitening® if I am wearing braces or if I have dental restorations (bridges, implants, crowns or veneers)?

In this case use the INTENSIVE mode to remove plaque and bacteria. In case of braces, place the brush gently between your brackets and gums and move the brush in a light circular motion along your teeth to keep both brackets and teeth free from debris.

I have recently had oral surgery.

Can I use Black Whitening®?

Always discuss your oral hygiene requirements with your oral surgeon or dentist and follow the given instructions after an operation. As soon as your teeth and gums heal to permit brushing, you can proceed using Sonic Generation®.

You will get further information at:

www.megasmile.com/sonic

BLACK WHITENING® BRUSH HEADS

**megasmile Sonic Black Whitening brush head
BW2 WHITENING Art. 200-223**

Swiss developed high-tech bristles with grains of activated charcoal neutralise dental plaque and polish teeth without scouring.

**megasmile SOFT Sonic Black Whitening brush head
BW2 SOFT Art. 200-224**

Swiss developed ultra-soft bristles with grains of activated charcoal neutralise dental plaque. The gentle massage effect is a relief for irritated gums.

CLEANING THE APPLIANCE

Rinse the brush head and bristles after each use with running water. Remove the brush head from the handle and rinse the connection at least once a week.

Always disconnect the charger from the socket, wipe it with a cloth and only plug it into the socket if fully dried. If you are not going to use the appliance for an extended period of time, remove the main plug from the wall socket. Then clean the appliance and store it in a cool and dry place.

WARRANTY

megasmile® guarantees its products for two years after the date of purchase. This warranty is the same in every country in which this device is sold by megasmile® or another authorized dealer.

This warranty does not cover:

Defects due to inappropriate use, regular wearing off or handling, especially of the exchangeable brush, but also flaws that insignificantly affect value or functionality of the device.

The warranty runs out if being repaired by unauthorized services or if opening the body housing.

In order to claim warranty, please send the entire device with proof of purchase to the official megasmile® partner in your country. If you do not know your official megasmile® partner, please visit our homepage

www.megasmile.com

or contact us directly in Switzerland:

+41 (0)44 520 0417

(The usual costs of a phone call into a foreign country will be charged).

TECHNICAL SPECIFICATIONS

Mains voltage: 100-240 V AC , 50 Hz, 2.5 W, Battery: rechargeable Li-Ion-battery

Operating time: 40 minutes (full to empty battery)

Initial charging time: 24 hours.

This device complies with the guidelines of the EU (EMC/RoHS) and low voltage.

Disposal

The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the power supply when removing the battery. Dispose of the discharged batteries correctly. Special containers are provided for disposing of

rechargeable batteries at shops selling them or public collection points. When disposing the main unit, affix scotch tape or label to the electrodes of the battery. Do not throw the rechargeable battery in fire.

Correct disposal of this product:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

megasmile® ist eine eingetragene Marken der Megasmile AG, Teufen, Schweiz.
Alle Rechte und Änderungen vorbehalten.

megasmile® est une marques enregistrée de la compagnie Megasmile AG,
Teufen, Suisse. Tous droits et changements reserves

megasmile® is an internationally registered trademarks of Megasmile AG,
Teufen, Switzerland. All rights and modifications reserved.

www.megasmile.com | hotline@megasmile.com

MS-BW4 Art. 200-220 SONIC BLACK WHITENING®