
Description	1027.2 LG0779.01 - Slim Sonic Leaflet - IEUK		
Date	31st July 2014	Proof N°	04
Perigord N°	1027.2	Component	Leaflet

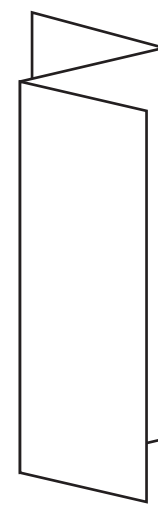3 PANEL
HARMONICA
FOLD

ENG Please read these instructions carefully before use

Congratulations on the purchase of your new **SLIM Sonic** toothbrush! **SLIM Sonic** operates at 22,000 brushstrokes per minute.

How to Use

- Put toothpaste on the brush.
- Press the power button firmly until vibrations begin.
- Brush teeth as recommended by your dentist.
- To turn **SLIM Sonic** off, press button firmly until vibrations stop.

How to Care for Your SLIM Sonic

- After every use, rinse the brush head thoroughly.
- Shake excess water from bristles and dry your brush head before covering and storing.

Avoid getting water on button.

- Remove the battery if you won't be using the **SLIM Sonic** for an extended period of time.
- To clean your **SLIM Sonic** toothbrush, wipe surface with a damp cloth. Make sure the battery box is secure to prevent water from getting into the motor or battery compartment.

Do not immerse your SLIM Sonic in water. Each SLIM Sonic brush head will last approximately 10 - 12 weeks.

Ordering extra brush heads:
Replacement **SLIM Sonic** brush heads are available to purchase.

Contact Details:

Lifes2good 7 Racecourse Business Park, Ballybrit, Galway City, Ireland Tel: 1890 20 26 20 www.lifes2good.com	Freepost RRAG-BAKK-JAYE Lifes2good UK Radlett, WD7 7AL United Kingdom Tel: 0845 601 8840
---	---

Limited 2 Year Warranty
The Violife LLC warranty lasts two years from the date of purchase.
Keep original receipt for warranty service. During this warranty period, return the product to your vendor, with proof of purchase.

Ireland:
Lifes2good,
7 Racecourse Business Park,
Ballybrit, Galway City, Ireland
Tel. 1890 20 26 20

UK:
Freepost RRAG-BAKK-JAYE
Lifes2good (UK) Ltd, Radlett, WD7 7AL, United Kingdom
Tel. 0845 601 8840

Important: This warranty does not cover damages resulting from accidents, misuse, alterations, abuse, lack of reasonable care, or use other than that described in these instructions, or performance of services or repairs by an unauthorised service centre.

WARNING:
This unit contains small pieces that may pose a choking hazard. Keep out of reach of children.

Changing the Replaceable Toothbrush Head (See pg.4)

- Firmly hold the **SLIM Sonic** base section with one hand and the base of the brush section with the other.

NOTE: If you do not rotate the brush head you risk damaging the brush head or the main body.

- Turn the brush section counter-clockwise and when the brush disengages, pull the head off the base.
- Place the new brush head on the stem of the base and turn clockwise until the brush is secure on the base.

Replacing the Battery (See pg.5)

Alkaline batteries are recommended

- Grip the **SLIM Sonic** at the base of the brush section with one hand and the bottom of the base section with the other.
- Turn the base section counter-clockwise and pull down, exposing the battery box.
- Remove the old battery. Insert a new AAA battery with the (-) end facing up.
- Replace the base cover and turn clockwise until the base section is securely in place.

NOTE: When battery is low, vibration gets weaker.

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A garantia Violife LLC. Tem a duração de dois anos a contar da data de compra.