PELVIC PARTNER
Security Notice Booklet

Technical Drawing 1.0

Ramblin Brands Pte Ltd
November 2025


Dimension: H 105 mm X W 99mm (Open)
H 105 mm X W 49.5 mm (Closed)

Construction: Facing pages booklet

Material: Uncoated white paper stock,
70 - 80 gsm (normal paper)

Finish: Matte

Colour: Black and white

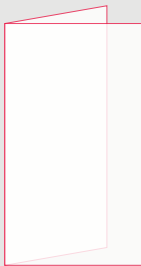smile😊makers

**INTIMATE ESSENTIALS
FOR HEALTHY HABITS**

Smile Makers makes self-care easy,
bringing comfort, confidence, and
connection to every part of your
routine – today and for the future.

English

EN

USAGE INSTRUCTIONS
 We are delighted you have chosen to try Pelvic Partner by Smile Makers. We hope you find the experience pleasant and pampering. We've tried to make it easy to use, but here is some guidance to help you along.

A light will blink when you need to charge your product. Once plugged in for charging, the light will remain lit until the battery is 100% charged. Always charge your product using only original chargers provided by us with your purchase. Before every charge, make sure the DC socket is free from water or lubricant residue.

Pelvic Partner has 1 button.

- To turn on the product, press the button with a single click.
- To turn the product off, hold the button for 2 seconds.
- To turn on/off the travel lock, press again for 4 seconds.
- Pelvic Partner has 8 progressive levels.
- Clicking the button takes you through the different modes.

This product is for use only by adults. If at any point you feel discomfort, nauseous, or otherwise unwell, please stop using the product immediately. If you are pregnant (congratulations!), have a pacemaker, diabetes, phlebitis or thrombosis, please consult with a medical professional regarding the risk of blood clots before using any personal massager.

CLEANING, STORAGE AND SAFETY
 The surface of Pelvic Partner is made of cryogenically-polished, FDA-approved food-grade silicone. So it is super safe, doesn't attract much dust or lint, and is really easy to clean.

Always clean Pelvic Partner thoroughly before and after each use. Wash with lukewarm water and pat dry with a lint-free cloth or towel. Never use cleansers containing alcohol, petrol or acetone. It should go without saying, but for insurance purposes we need to so here it is - never place your product directly in a pot of boiling water.

Pelvic Partner is waterproof, so it is safe to wash under running water or even in the shower. However, avoid fully submerging and getting water into the charging socket.

Never use Pelvic Partner together with silicone-based lubricants, massage oils or hand creams.

Avoid leaving any product from Smile Makers in direct sunlight and do not expose them to extreme heat. Do not use on swollen or damaged skin. Keep out of the reach of children.

In case of malfunction, please contact us at customercare@ramblinbrands.com and we will help you out.

WARRANTY TERMS & CONDITIONS
 Pelvic Partner comes with a warranty for a period of TWO (2) YEARS, after the date of original purchase, against defects due to faulty workmanship or materials. If you discover a defect and notify us during the warranty period, we will replace it free of charge.

The warranty covers working parts that affect the product's function. It does not cover cosmetic deterioration caused by fair wear and tear or damage caused by accident, misuse or neglect. Any attempt to open or take apart the product will void the warranty.

Claims under warranty must be supported by reasonable evidence that the date of the claim is within the warranty period. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period.

To claim your warranty, please email customercare@ramblinbrands.com informing us of the nature of the fault, and attaching evidence of purchase. In certain circumstances, you may be asked to return the product to us for inspection. This is to ensure we investigate faulty products and improve our processes for future users.

DISPOSAL INFORMATION
 Help us take care of the environment and do not dispose this product along with regular household waste. Please bring it to a collection point for recycling of electrical and electronic equipment.

PRODUCT SPECIFICATIONS
 IP Rating: IPX6
 Battery: 3.7V DC lithium
 Charging: USB2.0
 Max Voltage: 5.0V DC
 Frequency: 145±10% Hz
 Max noise level: ≤55 dB
 Interface: Level 1/2/3/4/5/6/7/8
 Weight: 78 g

Disclaimer: Ramblin' Brands reserves the right to revise this publication and to make changes from time to time in the contents hereof without obligation to notify any person of such revision or changes. The product exterior and technical design may be changed for improvements without notice.

© 2025 ALL RIGHTS RESERVED.
 Ramblin' Brands Pte. Ltd., 38 Beach Road #23-11, South Beach Tower, Singapore 189767.

Website: www.SmileMakersCollection.com
 Support: customercare@ramblinbrands.com

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©2025 ALL RIGHTS RESERVED.
Ramblin’ Brands Pte. Ltd., 38 Beach Road